

OPERATION AND SERVICE MANUAL FOR THE REFRACTIVE INDEX DETECTOR **RI2000**



SRI Instruments Europe GmbH
Menzenberger Str. 93 · D – 53604 Bad Honnef, Germany
Telephone: +49 (0) 2224 766 15 · Telefax: +49 (0) 2224 988 09 43
<https://sri-instruments-europe.com>

Version 1.0
Printed: 20.09.2023

File: *RI2000 english.doc*
Last modified: *KHS*

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1 CONTACT ADDRESS	4
2 PRECAUTION AND WARNINGS	5
3 ELECTRICAL WARNINGS	5
4 GENERAL WARNINGS	6
5 SPECIFICATIONS	8
6 PHYSICAL SPECIFICATIONS	9
7 GUARANTEE CONDITIONS	10
8 PARTS INCLUDED IN THE DELIVERY	11
9 FRONT PANEL AND KEYBOARD	12
9.1 THE DISPLAY	12
9.2 THE AUTOZERO BUTTON	13
9.3 THE PURGE BUTTON	13
9.4 THE BUTTON POL +/-	13
9.5 THE BUTTONS ARROWLEFT (◀) AND ARROWRIGHT (▶)	13
9.6 THE BUTTONS ARROWUP (▲) AND ARROWDOWN (▼)	13
9.7 THE MENU BUTTON	14
9.8 THE ENTER BUTTON	14
9.9 THE IN PORT	14
9.10 THE OUT PORT	14
10 THE BACK SIDE OF THE RI2000	15
10.1 OUTPUT-CONNECTORS	15
10.1.1 THE INTEGRATOR SIGNAL OUTPUT	15
10.1.2 THE RECORDER SIGNAL OUTPUT	15
10.1.3 THE DIGITAL OUTPUT	15
10.2 INPUT-CONNECTORS	16
10.2.1 THE EXTERNAL START SIGNAL	16
10.2.2 THE EXTERNAL AUTOZERO SIGNAL	16
10.2.3 THE EXTERNAL PURGE SIGNAL	16
10.3 THE RS232 PORT	16
10.4 MASS SCREW	17
10.5 OPTICAL ADJUST	17
10.6 DRAIN	17

10.7	MAINS ADAPTOR	17
10.8	FUSES	17
11	<u>INSTALLATION OF THE RI2000</u>	18
11.1	CONNECTING THE DETECTOR TO THE CHROMATOGRAPHY SYSTEM	18
11.2	CONNECTION TO AN ANALOG DATA SYSTEM	19
11.3	CONNECTION TO A DIGITAL DATA SYSTEM (PC)	19
12	<u>PRINCIPLE OF DETECTION</u>	20
12.1	THE OPTICAL SYSTEM OF THE RI2000	20
13	<u>FLOW PATHS IN THE RI2000</u>	21
13.1	MEASURING MODE	21
13.2	PURGE MODE	22
14	<u>OPERATING THE RI2000 DETECTOR</u>	24
14.1	OVERVIEW FIRMWARE VERSION 3.2C	24
14.2	THE INITIALIZATION STAGE	25
14.3	THE NORMAL MEASURING MODE (NORMMODE)	26
14.4	TEMPERATURE SETTINGS	28
14.5	SETTING SIGNAL SMOOTHING	30
14.6	SETTING RECORDER RANGE – SIGNAL AMPLIFICATION	31
14.7	SETTING RECORDER OFFSET AND RECORDER EXTENT	32
14.7.1	RECORDER OFFSET	32
14.7.2	RECORDER EXTENT	33
14.8	MARKERS	34
14.9	SETTINGS FOR SERIAL PORT (RS232)	35
14.10	MENU STRUCTURE IN NORMAL MEASURING MODE (NORMMODE)	36
15	<u>THE SERVICE MODE</u>	37
15.1	FIRMWARE INFORMATION	38
15.2	VIEWFINE MODE	39
15.3	ADJUSTING DIGITAL – ANALOG – CONVERTERS	41
15.4	CALIBRATING THE RI2000	42
15.4.1	HOW TO CALIBRATE...	43
15.5	SERVMODE OVERVIEW	46
15.6	SERIAL COMMUNICATION	48
16	<u>MAINTENANCE OF THE RI2000</u>	49
16.1	IDENTIFYING HARD- AND SOFTWARE	49
16.2	ERROR MESSAGES	50
16.2.1	MALFUNCTION OF THE TEMPERATURE SENSOR	50
16.2.2	MALFUNCTION OF THE ANALOG DIGITAL CONVERTER	51
16.2.3	INTENSITY ALARM	51

16.2.4	SIGNAL OUT OF RANGE	52
16.2.5	RESET CONFIGURATION	53
16.3	ADJUSTING THE OPTICAL SYSTEM	55
16.4	LAMP EXCHANGE AND ADJUSTMENT OF THE LAMP	57
16.5	CHECKING SUM- AND DIFFERENCE VOLTAGES	61
16.6	CHECK AND REPLACEMENT OF THE VALVE	64
16.7	CHECKING AND CLEANING THE FLOW CELL	66
16.8	THE HEATING CIRCUIT OF THE RI2000	69
17	<u>VOLTAGE TEST POINTS</u>	70
18	<u>SPARE PARTS</u>	71
19	<u>TEST POINTS FOR LAMP VOLTAGE (FOR SERIAL NUMBER 0312XXX)</u>	72
20	<u>TEST POINTS FOR LAMP VOLTAGE (FOR SERIAL NUMBER 0401XXX)</u>	73
21	<u>VOLTAGE TEST POINTS (OVERVIEW) FOR 0312XXX SERIES</u>	74
22	<u>VOLTAGE TEST POINTS (OVERVIEW) FOR 0401XXX SERIES</u>	75
	<u>DECLARATION OF CONFORMITY</u>	76

1 CONTACT ADDRESS

In case of problems with your RI2000 refractive index detector please contact your local distributor or SRI Instruments Europe GmbH.

SRI Instruments Europe GmbH
Menzenberger Str. 93
D – 53604 Bad Honnef
Germany

Phone: +49 (0) 2224 766 15
Fax: +49 (0) 2224 988 09 43

<https://sri-instruments-europe.com>
mail: mailto@sri-instruments-europe.com



[GC](#) [HPLC](#) [Downloads](#) [Company](#) [Contact](#)

Call us: +49 (0)2224 766 15



Welcome

SRI Instruments Europe GmbH is a worldwide manufacturer and distributor of gas chromatographs, autosamplers, HPLC systems, software, columns and accessories.

Our scientific partners are institutes, universities and authorities as well as laboratories from the environmental – food – and biotech industry since 30 years.

Innovative electronics and a good reproducibility of the analysis results are our trademark. Our customers can expect tailor-made solutions and service-oriented action.

We are characterized by a high reliability – in case of problems we will not leave you alone.

Product Highlights



[GC Catalog](#)



[HPLC Catalog](#)

2 PRECAUTION AND WARNINGS

The manufacturer is not reliable for any damage, harm or financial loss caused directly or indirectly by the use of this instrument if the instrument is handled without the observation of this manual or handling with carelessness.

■ **NOTE**

Please read this manual carefully before working with the RI2000 refractive index detector. In case of any question please do not hesitate to contact you local distributor or SRI Instruments Europe GmbH directly.

3 ELECTRICAL WARNINGS

- Before opening the housing of the instrument make sure the detector is switched off and disconnected from power supply.
- The voltage selected at the fuse on the backside of the instrument has to be set correctly to 110 V or 220 V according to your local power supply.
- Before changing the selected voltage the power supply has to be disconnected.
- The detector may only be connected to plugs with grounding.
- The RI2000 may only be operated in connection with other devices which fit to the safety requirements.

4 GENERAL WARNINGS

- To prevent damages of the RI2000 all capillaries and cables have to be checked for damages and leakages.
- For the disposal of inflammable and/or toxic solvents a plan for waste management has to be created. Such solvents may not get into the drain.
- The RI2000 refractive index detector is built to operate at temperatures between 10 °C and 35 °C.
- To reach a reliable operation of the detector it is recommended to use filtered samples and solvents only.
- The detector may be cleaned with appropriate cleaning agents only.
- Make sure that no liquid gets inside the detector. Liquid inside the housing may cause electrical short circuits which may result in the damage of the instrument.
- To prevent electrical shocks make sure that the detector is disconnected from power supply when the housing is opened to perform service work inside.
- To disconnect the instrument from power supply simply unplug the power cable.
- Electronic circuit boards and electronic components are sensitive to electrostatical charges.
- For some maintenance operation it is required to open the instrument's housing. Make sure that the instrument is disconnected from power supply before removing the housing. If it is necessary to remove the housing during operation make sure not to touch electrical parts inside the detector.

- In case of the use of dangerous solvents pay attention to safety instructions regarding this solvent.
- Solvents should be degassed before use with a refractive index detector.
- After use of salt containing solvents (such as buffers) the detector should be purged with distilled water.
- Make sure that the 6 bar (0.6 MPa) pressure limit of the flow cell is not exceeded.
- Make sure that the 2 bar (0.2 MPa) pressure limit of the valve is not exceeded.
- During operation the housing of the detector should be closed.
- Do not use the detector in ambience of aggressive gases, very high humidity, strong vibrations and strong changes in the ambience's temperature.

5 SPECIFICATIONS

Type of detection:	Refractive Index
Effective range for refractive index:	1.00 – 1.75 RIU
Effective range:	$\pm 1000 \mu\text{RIU}$
Optical null balance:	In the whole effective range by adjusting the mirror
Signal null balance:	In the whole effective range by using the AutoZero function
Materials which contact the mobile phase:	Stainless steel, glass, PTFE
Maximum pH-range:	2.3 – 9.5
Temperature control for optical bench:	7°C above ambient temperature, 35 °C up to 55 °C in 1 °C steps
Analog signal output:	Recorder / integrator
Noise (analog signal):	$\pm 7 \text{ nRIU}$ ($\pm 7 \mu\text{V}$)
Noise (digital output):	$\pm 3 \text{ nRIU}$ ($\pm 3 \mu\text{V}$)
Drift:	$< 1 \mu\text{RIU/h}$ ($< 1\text{mV/h}$)
Integrator output:	$\pm 1.2 \text{ V}$ (fix)
Recorder output:	$\pm 1.2 \text{ V}$ (adjustable)
Recorder range:	8 steps selectable in the region from 1/8 to 16/1
Recorder offset:	Selectable: 0 mV, 10 mV, 100 mV
Signal range (recorder):	Selectable: 10 mV, 100 mV, 1000 mV
Recorder marker:	Marker function ON/OFF
Flow cell:	Glass cell (quartz) with two chambers
Volume:	9 μL
Angel:	45 °
Max. Pressure:	6 bar (0.6 MPa)
Max. Flow rate:	0.1 – 3.0 mL/min
Valve:	3/2 path valve, 12 V
Max. Pressure:	2 bar (0.2 MPa)
Volumes:	Sample entry → sample cell: 24 μL Sample entry → sample cell → reference cell → sample exit: 1500 μL
Communication:	Digital interface: RS232, bidirectional Digital input: TTL (Purge, AutoZero, Start/Marker) Digital output: TTL (Intensity alarm)

6 PHYSICAL SPECIFICATIONS

Dimensions: 220 mm * 165 mm * 350 mm (W * H * D)

Weight: 12 kg

Voltages: Mains adaptor with voltage selector
100 – 120 V
220 – 240 V

Frequency: 50 or 60 Hz

Power consumption: max. 50 W

Temperature range: 10 °C to 50 °C

7 GUARANTEE CONDITIONS

The term of guarantee depends on your local law. Besides this SRI Instruments Europe GmbH affords guarantee for at least 12 months beginning after purchase of the RI2000 refractive index detector.

All instruments are tested and certified by SRI Instruments Europe GmbH quality control.

Only defects which result from faulty manufacture or material defects are covered by the guarantee. In case of a defect the original sales slip will be needed to make use of the guarantee. Repair works covered by the guarantee may only be performed by SRI Instruments Europe GmbH or licensed distributors.


The following cases are not covered by the guarantee:

- Improper use (e. g. capacity overload, use of not approved tools) of the RI2000 detector
- Parts which are liable to aging or abrasion like lamps, valves, heater cartridges or flow cells
- Damages caused by use of force or not approved tools
- Damages which result from improper use due to non-observance of the operation manual, the use under abnormal conditions or improper maintenance
- Defects resulting from the use of third-party parts which are not approved by SRI Instruments Europe GmbH
- Instruments which are modified after purchasing by the customer
- Normal abrasion
- Fully or partial disassembled refractive index detectors

8 PARTS INCLUDED IN THE DELIVERY

Quantity	Description
1 EA	Refractive Index Detector RI2000
1 EA	Stainless steel capillary, 0.25 mm inner diameter, 1.59 mm outer diameter (for sample entry)
1 EA	Stainless steel capillary, 1.00 mm inner diameter, 1.59 mm outer diameter (for sample exit)
1 EA	PTFE tube, 1.59 mm inner diameter, 2.80 mm outer diameter (for sample exit)
2 EA	Screw connection (for sample entry and sample exit)
1 EA	Power cable
2 EA	Fuses (0.5 A, slow)
2 m	Signal cable, 2-wired
1 EA	Hex-wrench (for adjustment of the mirror, \varnothing 2.5 mm)
1 EA	Operation & Service Manual

9 FRONT PANEL AND KEYBOARD

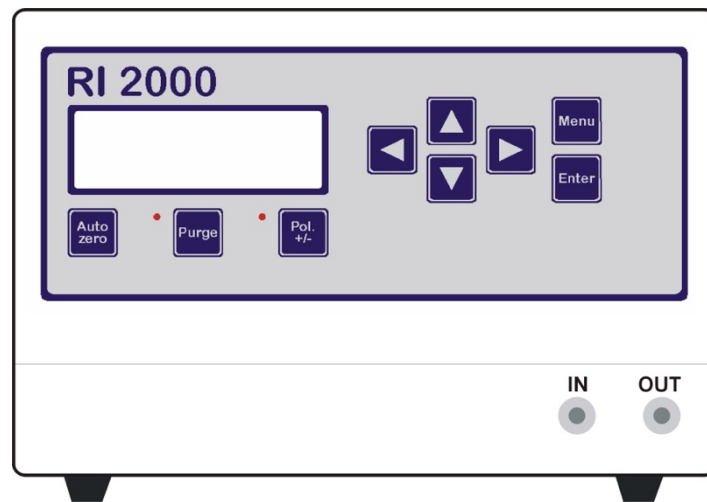

9.1 THE DISPLAY

The RI2000 refractive index detector is equipped with a liquid crystal display to show the current detector signal as well as system parameters. This display contains four lines with 20 characters each.

```

Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
NormMode OK Menu ,◇,↙
  
```

Line 1 shows the current temperature of the optical bench (in °C).

Line 2 shows the optical balance (in %%).

Line 3 shows the current detector signal (in mV).

The value shown here is identical to the digital output signal and to the analog signal which can be recorded at the recorder output.

Line 4 This line shows the status of the RI2000 detector. Every few seconds the displayed parameter changes. By this you can get a complete overview of all parameters in a short time.

9.2 THE AUTOZERO BUTTON

By pressing the **AutoZero** button the detector signal is internally set to Zero. To return to the formal (unmodified) signal, press the **AutoZero** button for approx. three seconds.

9.3 THE PURGE BUTTON

To flush the reference chamber of the flow cell with mobile phase press the **PURGE** button once. If the purge mode is activated a red LED will light up. Press the **PURGE** button once again to deactivate the purge mode. Now the mobile phase will flow only through the sample chamber of the flow cell.

9.4 THE BUTTON POL +/-

By pressing this button the polarity of the detector signal will be changed. A negative signal like e. g. -5mV will be transformed into a positive signal $+5\text{mV}$ when the polarity mode is activated. A red LED will light up when this mode is activated.

9.5 THE BUTTONS ARROWLEFT (◀) AND ARROWRIGHT (▶)

This arrow buttons have several functions. If only one button is pressed you can change the settings of one parameter to smaller (◀) or bigger (▶) values. If both buttons are pressed together for about three seconds the detector will be set to the service mode.

9.6 THE BUTTONS ARROWUP (▲) AND ARROWDOWN (▼)

If different options for one parameter are available you can select the desired option by using these arrow buttons. **ArrowUp** (▲) will switch to higher values and **ArrowDown** (▼) to smaller values.

If the marker function is activated you can send a marker signal to the recorder output by pressing the **ArrowUp** (▲) button (see chapter 14.8).

9.7 THE MENU BUTTON

By pressing this button you can select a special sub menu or discard changes you performed on a parameter.

9.8 THE ENTER BUTTON

Press the **ENTER** button to open the selected sub menu and to accept changes you performed on a parameter.

9.9 THE IN PORT

This is the entrance for the mobile phase coming from the column. Connect the capillary coming from the column to this port using screw connectors.


9.10 THE OUT PORT

This is the exit of the mobile phase. Connect you waste bottle to this port using the PTFE tube.

NOTE

Use only capillaries with an inner diameter of 1 mm for the exit of the mobile phase. By this blockage can be prevented. Make also sure that the pressure limits for the flow cell (6 bar, 0.6 MPa) and the valve (2 bar, 0.2 MPa) are not exceeded.

10 THE BACK SIDE OF THE RI2000

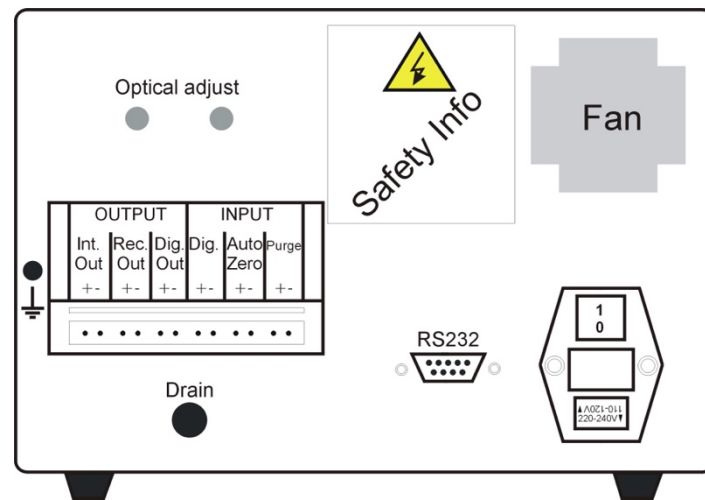

10.1 OUTPUT-CONNECTORS

The RI2000 refractive index detector has two analog and one digital signal output.

10.1.1 THE INTEGRATOR SIGNAL OUTPUT

If you want to record your data using an analog integrator device, connect the integrator using a shielded signal cable to the integrator output (Int. Out). Make sure that the polarity (+/-) is correct.

10.1.2 THE RECORDER SIGNAL OUTPUT

To record data using an analog recorder, connect the device using the shielded signal cable to the recorder output (Rec. Out). Make sure that the polarity (+/-) is correct.

10.1.3 THE DIGITAL OUTPUT

The RI2000 gives a TTL pulse in case of an intensity alarm. This pulse is given to the digital output (Dig. Out). If you want to record the digital detector signal you need to use the RS232 serial connector.

10.2 INPUT-CONNECTORS

The RI2000 refractive index detector has three digital (TTL) inputs. These can be used to control the system externally.

10.2.1 THE EXTERNAL START SIGNAL

Use the digital input connector (Dig.) to connect an external start signal (e. g. generated by an autosampler or a manual injector) to the instrument.

10.2.2 THE EXTERNAL AUTOZERO SIGNAL

Use this input connector (AutoZero) to connect an external AutoZero signal to the instrument.

10.2.3 THE EXTERNAL PURGE SIGNAL

Using this connector (Purge) you can use an external signal to activate the detector's purge mode.

NOTE

All TTL signals handled by the RI2000 detector use a TTL level which is based on instrument's ground potential.

10.3 THE RS232 PORT

The RS232 port can be used to acquire the detector's signal using a personal computer as well as to control the instrument using the computer. To communicate with the detector using a computer you need to use a certain protocol which is explained in [chapter 15.6](#) of this manual.

10.4 MASS SCREW

Use the mass screw to connect the shielding of the analog signal cable to the detector's chassis. This will result in a smoother detector signal.

10.5 OPTICAL ADJUST

The two black caps on the back side of the detector cover holes in the detector's housing. If needed the optical system of the detector can be adjusted using the hex-wrench which came with the instrument. To adjust the mirror, remove these caps and turn the adjustment screws *carefully*.

10.6 DRAIN

In case of a leakage inside the housing mobile phase can leave the housing using the tube. You might connect an additional waste tube to this exit.

10.7 MAINS ADAPTOR

Use the power cable which was delivered with the detector to connect the instrument to the power supply. Make sure that the correct voltage is selected before connecting to the mains. To high voltage will damage the instrument.

10.8 FUSES

Next to the mains plug you can find the instrument's fuses. The type of this fuse depends on the voltage the instrument is operated.

220 V:	2 x 500 mA, slow
110 V:	2 x 1000 mA, slow

Inside the housing is an additional internal fuse (1 x 500 mA, slow).

11 INSTALLATION OF THE RI2000

11.1 CONNECTING THE DETECTOR TO THE CHROMATOGRAPHY SYSTEM

Connect the INLET on the front of the detector to the capillary which comes from the column of your chromatography system. To attach the column to the detector use a 1/16" stainless steel capillary with 0.25 mm inner diameter.

Connect the 1/16" stainless steel capillary with 1 mm internal diameter to the PTFE tube and the OUTLET port of the instrument. Use the screw connections to fix the capillaries to the INLET and OUTLET port.

Place the free end of the PTFE tube in a bottle for solvent waste. The droplets leaving the tube should run along the bottle's wall and not fall down. The waste bottle and the refractive index detector should be positioned at the same high level.

If you need a solvent mixture to perform your analytical separation it is recommended to mix the different components manually before using the solvent. The use of a gradient pump to generate a solvent mixture is not possible as the mixture generated by the gradient pump is not consistent. Small changes in the solvent compositions will result in unstable detector signals.

It is also recommended to degass the solvents before use with a refractive index detector as this will result in a more stable detector signal.

■ **NOTE**

Before dispatch of the instrument the flow cell was purged first with ethanol followed by air. However, there are still some rests of ethanol in the flow cell. For transportation or long time storage of the instrument the procedure of purging with ethanol and air is recommended.

■ **NOTE**

The flow cell is built of specialized optical glass which is sensitive to pressure. The maximum pressure for the flow cell is **6 bar (0.6 MPa)**. If you want to use several detectors in series, make sure that the refractive detector is the last one in the row to prevent damages of the flow cell due to high pressure.

11.2 CONNECTION TO AN ANALOG DATA SYSTEM

Use a small screwdriver to connect the signal cable to one of the two pin connectors which are delivered with your instrument. Connect the signal cable either to the integrator output (Int. Out) or to the recorder output (Rec. Out) of the detector if you want to use an analog integrator or recorder as well as an analog data system to collect data. Make sure that the polarity (+/-) of the connection is correct.

■ NOTE

The signal coming from the integrator output (Int. Out) is a voltage between -1 V and 1 V. In case of a signal gain of 1 set in the detector parameters the signal at the recorder output (Rec. Out) is the same.

If you need a very high sensitivity it is recommended to use the recorder output (Rec. Out) as you can use the instrument's internal signal amplifier which might amplify the detector signal up to 16 times.

11.3 CONNECTION TO A DIGITAL DATA SYSTEM (PC)

The RS232 port on the back side of the RI2000 refractive index detector can be used for data acquisition as well as instrument control. Using the RS232 port for external instrument control you might activate / deactivate the purge mode, send an external start / marker signal and AutoZero signal.

To communicate with the RI2000 detector using the serial interface you will need a special software package which utilizes a special communication protocol. Details of this protocol you will find in [chapter 15.6](#) of this manual.

There is a software package available which allows the data acquisition with the RI2000 detector. The data can be saved in user defined file formats for later data treatment using chromatography software packages like *PeakSimple* or *APEX*.


12 PRINCIPLE OF DETECTION

The RI2000 refractive index detector acts as a differential refractometer which measures the deflection of a light beam as a result of different refractive indices of solutions in the reference and the sample chamber of the flow cell.

The beam of a tungsten lamp passes a convex lens and two slit apertures (called slit 1 and slit 2) before reaching a second convex lens. After this the light passes the flow cell which is diagonally separated into two chambers (sample and reference chamber). Behind the flow cell a mirror is placed which reflects the beam through the flow cell, the second convex lens and the second slit aperture (slit 2) toward a light sensor.

The light sensor is placed under the light source and consists of two photo diodes. The both diodes generate a current which is proportional to the light intensity reaching the sensor's surface.

12.1 THE OPTICAL SYSTEM OF THE RI2000

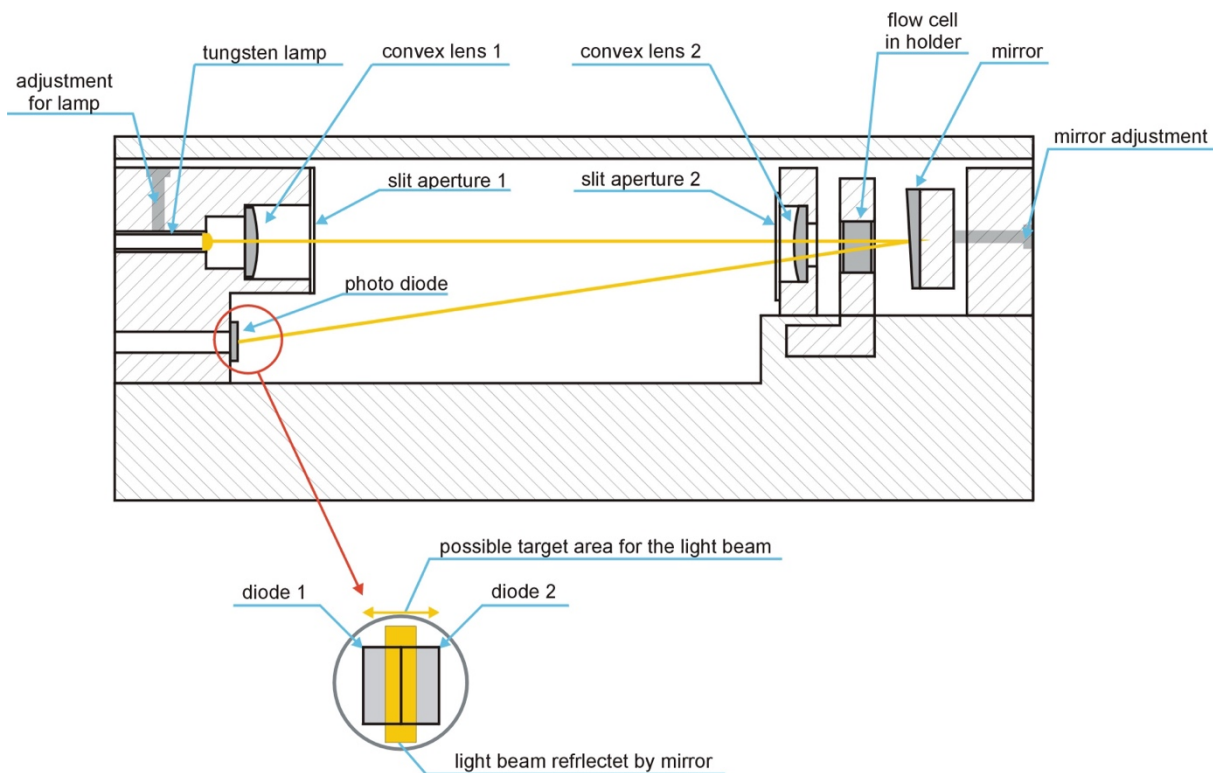


13 FLOW PATHS IN THE RI2000

13.1 MEASURING MODE


If the purge mode is disabled (red LED is off) the valve is switched that way that the mobile phase does not pass the reference chamber of the flow cell. The mobile phase comes from the column and enters the optical system through the INLET port and a heat exchanger to ensure that the measurement is performed at constant temperatures as the refractive index is highly dependent on the temperature. The temperature for the optical bench can be selected in the range from 35 °C to 55 °C. After passing the heat exchanger the mobile phase reaches the sample chamber of the flow cell. Leaving the flow cell the mobile phase flow again passes the heat exchanger going towards a T-connector which is connected with the purge valve. After passing the valve the liquid leaves the detector passing the OUTLET port.

The following figure illustrated the flow path in measuring mode with MTV - Valve.

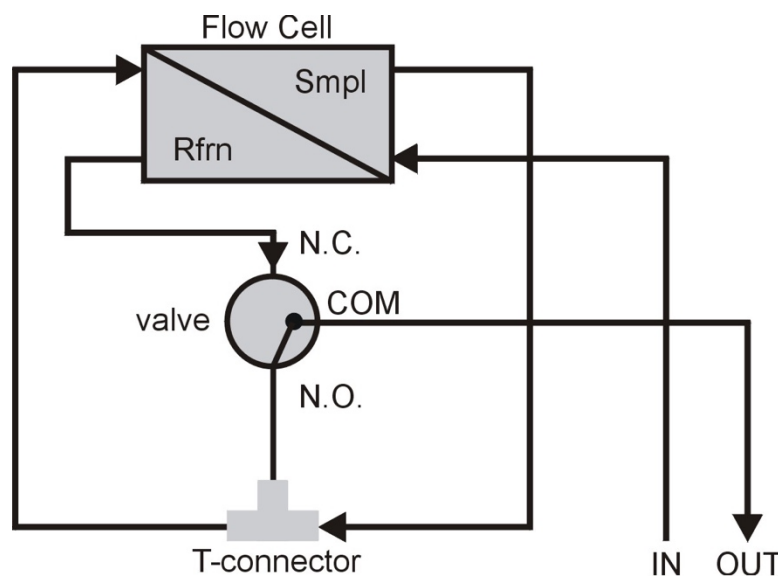

The following figure illustrated the flow path in measuring mode with ASCO - Valve.

For ASCO – valve the RI2000 - Firmware 4.2 is essential.

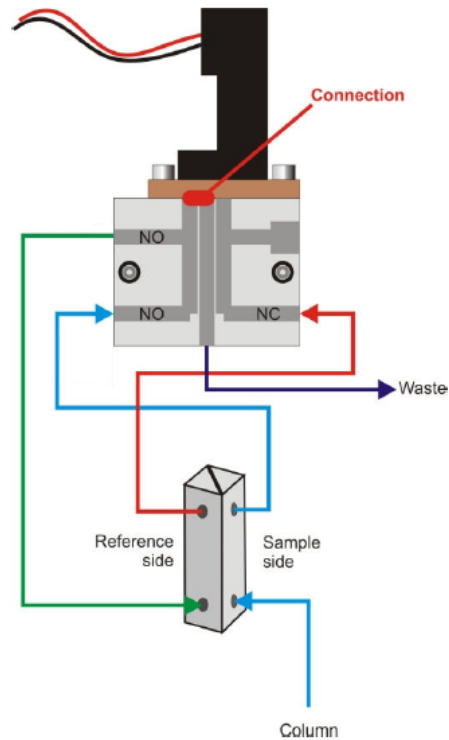

13.2 PURGE MODE


The purge mode is used to flush the reference chamber or the flow cell with fresh mobile phase (solvent). If the purge mode is activated the red LED will light up. The mobile phase now passes the reference chamber of the flow cell.

After leaving the column the phase enters the detector passing the INLET port. After this the mobile phase enters the optical bench passing the heat exchanger going towards the flow cell. Now the stream enters the sample chamber of the flow cell. Leaving the flow cell, passing the heat exchanger again the liquid reaches the T-connector and after this it reaches the flow cell. Now the reference chamber is flushed by mobile phase. Leaving the flow cell the liquid is guided through the OUTLET port out of the detector.


The purge mode should be activated regularly to ensure that the reference chamber contains mobile phase (solvent) which is similar to the mobile phase passing the sample chamber.

After purging the detector for a certain time you can start your measurement after deactivation of the purge mode by pressing the **PURGE** button once again.

The following figure illustrates the liquid stream in purge mode with MTV - valve.

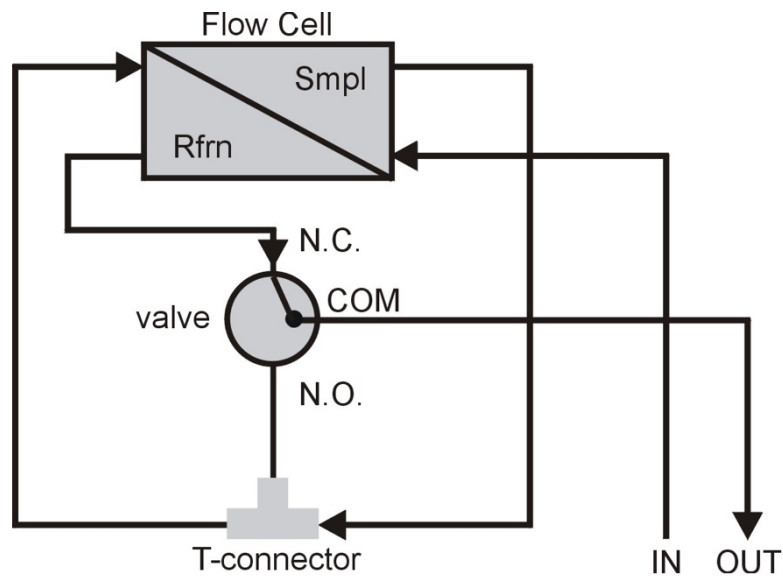

The following figure illustrates the liquid stream in purge mode with ASCO - valve.
For the new ASCO – valve the RI2000 - Firmware 4.2 is essential.

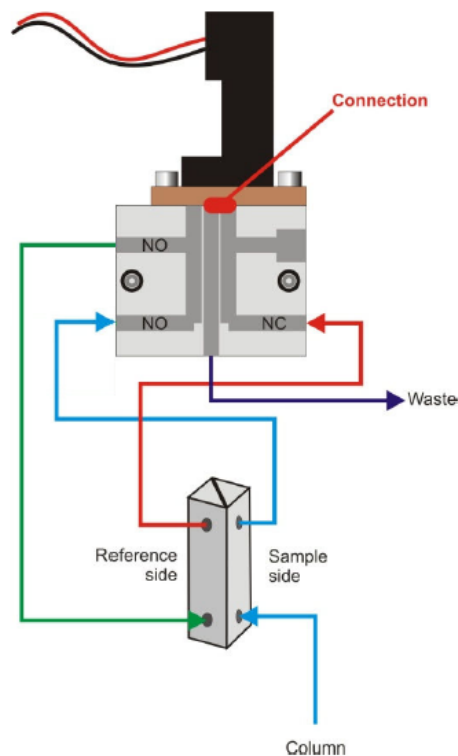


The purge mode should be activated until a stable baseline is reached.

After switching back to the measuring mode it might be necessary to wait once more until the baseline becomes stable.

14 OPERATING THE RI2000 DETECTOR

14.1 OVERVIEW FIRMWARE VERSION 3.2C

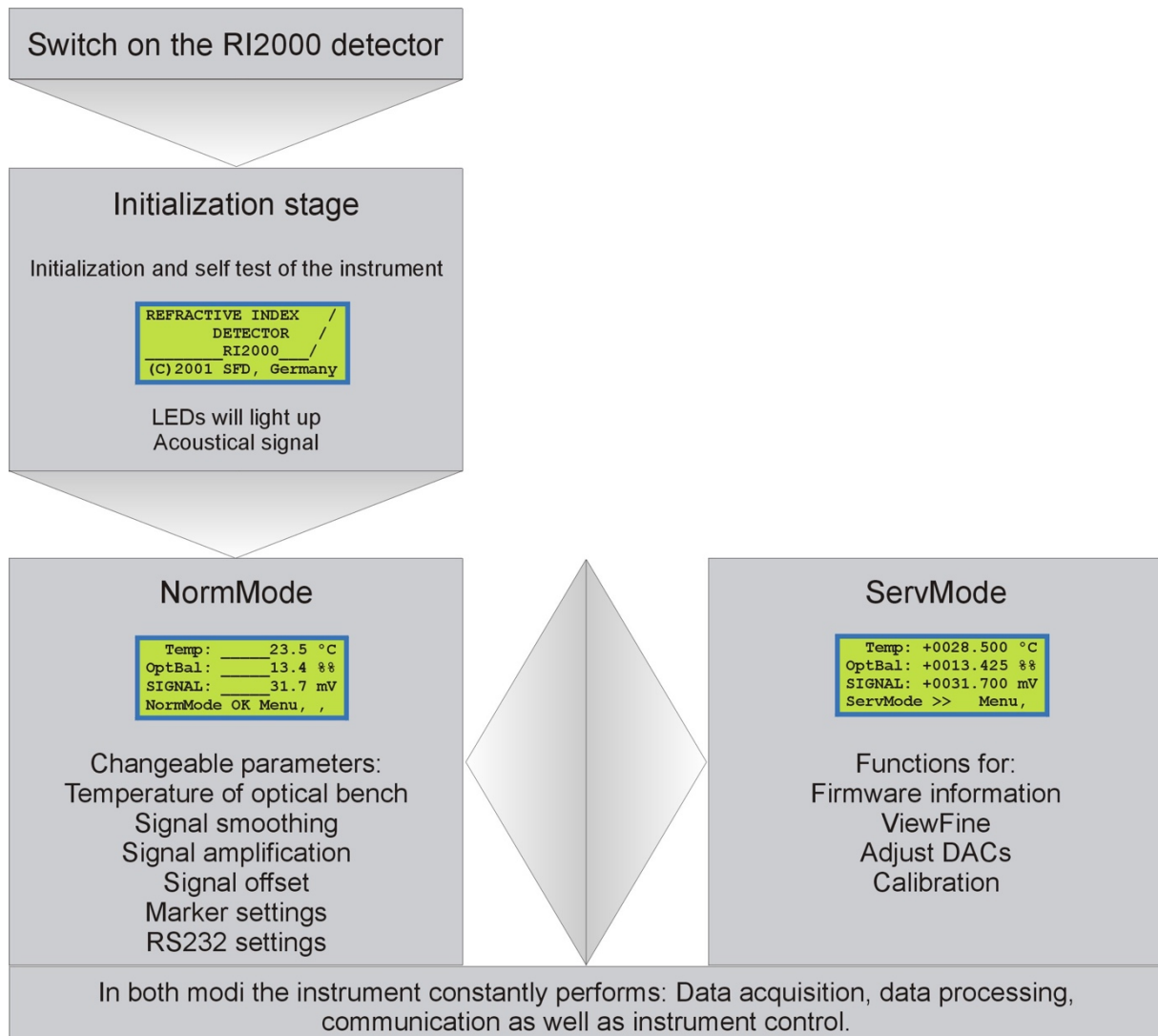

14.2 THE INITIALIZATION STAGE

The RI2000 refractive index detector does not consume any energy when it is switched off. The mains switch can be found on the back side of the instrument. After switching on the detector starts to initialize the internal micro controller and a testing procedure starts to check several system parameters. During this test the logo is shown in the display, the LEDs will light up for a short time and an acoustical signal will occur. The initialization procedure will be finished after approx. six seconds and the detector will switch automatically to the normal measuring mode (NormMode).

If you want to check the functionality of the internal keypad you can perform two different tests while the logo is shown:


1. Shortly press one of the buttons on the keypad. The detector will switch to normal measuring mode (NormMode) without waiting.
2. Press and hold one of the buttons. In this case the detector will not proceed to normal measuring mode. Instead of this it will show the logo until another button is pressed.

14.3 THE NORMAL MEASURING MODE (NORMMODE)

After performing the internal test procedure the detector will switch automatically to the normal measuring mode (NormMode). Now internally several tasks are running:

- Data acquisition
- Data processing
- Data output
- Internal instrument control

The display will show the following information:

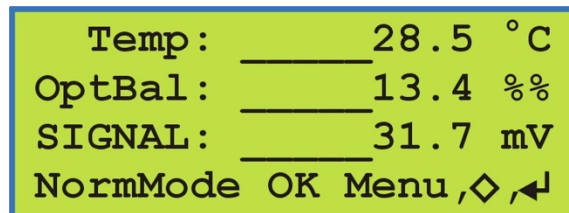

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
NormMode OK Menu ,◇,↵
```

Line 1 shows the current temperature of the optical bench in [°C].

Line 2 shows the optical balance between the two photo diodes of the light sensor in [%%] which corresponds with the relative position of the light beam on the surface of the light sensor.

Line 3 shows the current detector signal or error messages (in case of intensity alarms). The shown signal corresponds with the signal which can be measured at the integrator output (IntOut). In the normal measuring mode (NormMode) the signal is displayed with one decimal.

Line 4 is the status line and shows the current status of different detector parameters. Several parameters are shown. Every few seconds the display will change that you will have a complete overview of all parameters in a short time.

The following parameters are shown:

- SetHeating, °C Temperature control
- SigSMOOTH Signal smoothing, internally
- RecRANGE Recorder output settings
- RecOFFSET, mV Signal offset for recorder output
- RecMARKER Marker settings
- ComRS232 Settings for serial RS232 port

Another function of the status line is the input of new parameter settings by the user. The symbol ❖ represents the four ARROW keys (◀▶▼▲), ↵ stands for the **ENTER** button.

Using the menu button you can choose between several sub-menus which are accessed when pressing the **ENTER** key after selection. To adjust the settings for the desired parameters use the ARROW buttons. To accept new settings press the **ENTER** button. To discard new settings and keep the old value press the **MENU** button.

14.4 TEMPERATURE SETTINGS

Since the refractive index is highly depending on the temperature the temperature of the optical bench of the RI2000 refractive index detector can be controlled by the micro controller. To set the desired temperature press the **MENU** button until the following status line is displayed.

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
SetHeating, °C      OFF
```

Currently the temperature control is switched OFF. The optical bench is operated at ambient temperature plus additional 7 °C caused by internal heating resulting from different electronic components which heat up when operated. As a result of this the temperature can be controlled in a range from 35 °C up to 55 °C only. To change the temperature setpoint press the **ArrowUp** (▲) to switch to higher temperatures or use the **ArrowDown** button (▼) to select a lower temperature. If you set the temperature control to OFF the optical bench is operated at ambient temperature plus 7 °C.

Pressing the **ENTER** button accepts the new setting. In the normal measuring mode an asterisk (*) in the left top corner of the display indicates the activity of temperature control.

```
* Temp : _____ 35.0 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
SetHeating, °C      +35
```

The temperature sensor is constantly controlled by the internal firmware of the instrument. It is not possible to activate the heating when one of the following reasons occurs:


- The temperature sensor does not work or is not connected. In this case the error message “noTS” will appear in the display.
- The current temperature is below 9 °C or above 64 °C. In this case the error message “!” will be shown in the display and a periodic acoustic signal occurs. The acoustic signal will stop when the temperature reaches the range between 9 °C and 64 °C.

■ NOTE

Should the temperature control fail due to an electronic malfunction the heating will be switched off automatically at 72 °C by a thermo fuse ([see chapter 16.8](#)).

14.5 SETTING SIGNAL SMOOTHING

To adjust the grade of smoothing of the detector signal press the menu button until the following status line is shown in the display.

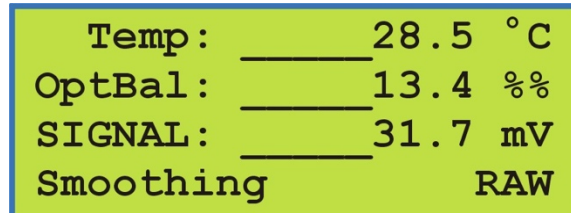

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
Smoothing                RAW
```

There are four different modes for signal smoothing:

- **RAW** No signal smoothing
- **FAST** Signal smoothing by calculating the average in a time interval of 0.4 seconds
- **MEDIUM** Signal smoothing by calculating the average in a time interval of 0.8 seconds
- **SLOW** Signal smoothing by calculating the average on a time interval of 1.6 seconds

The **ArrowUp** button (▲) changes to the next faster way of signal smoothing while the **ArrowDown** button (▼) switches to the next slower mode. To accept the new settings press the **ENTER** button.

14.6 SETTING RECORDER RANGE – SIGNAL AMPLIFICATION

To adjust the signal amplification press the **MENU** button until you can see the following status line in the display.

Temp :	_____	28.5 °C
OptBal :	_____	13.4 %%
SIGNAL :	_____	31.7 mV
RecRANGE		1/1

If the mode RecRANGE is activated it allows to amplify or attenuate the detector signal at the recorder output (RecOut).

If the amplification mode is set to “1/1” the signal at the recorder output (Rec. Out) corresponds to the signal at the integrator output (Int. Out).

By pressing the buttons **ArrowUp** (**▲**) you can switch to the next stronger amplification mode of the signal:

$1/1$ **▲** $2/1$ **▲** $4/1$ **▲** $8/1$ **▲** $16/1$

Using the **ArrowDown** (**▼**) button you can change the amplification to lower grades:

$1/1$ **▼** $1/2$ **▼** $1/4$ **▼** $1/8$

To accept your settings press the **ENTER** button.

14.7 SETTING RECORDER OFFSET AND RECORDER EXTENT

In this sub menu you can adjust the recorder offset or the recorder extent. If the recorder adjust is activated a certain (defined) voltage will be added to each data point recorded by the detector. By setting the recorder extent you may define a maximum for the detector signal.

To adjust settings press the **MENU** button until the following status line is shown in the display.

Temp :	_____	28.5 °C
OptBal :	_____	13.4 %%
SIGNAL :	_____	31.7 mV
OutREC ,mV		0add

By pressing the buttons **ArrowUp** (▲) and **ArrowDown** (▼) you can either add a static voltage to each data point recorded by the detector (detector offset) or define a maximum voltage for the detector signal (detector extent).

14.7.1 RECORDER OFFSET

You can define a static voltage of 0 mV, 10 mV or 100 mV which is added automatically to the displayed detector signal. This is reasonable in the case of negative baseline drift in connection with a data system which cannot handle negative detector signals. To adjust the recorder offset press the buttons **ArrowUp** (▲) and **ArrowDown** (▼) until the desired voltage is shown in the status line of the display.

Temp :	_____	28.5 °C
OptBal :	_____	13.4 %%
SIGNAL :	_____	31.7 mV
OutREC ,mV		100add

0 mV ▲ 10 mV ▲ 100 mV or 100 mV ▼ 10 mV ▼ 0 mV

To accept the new setting, press the **ENTER** button. To leave the mode without accepting press the **ENTER** button.

14.7.2 RECORDER EXTENT

To adjust the recorder extent press the **ArrowUp** button (▲) until the following status line is displayed.

Temp :	_____	28.5 °C
OptBal :	_____	13.4 %%
SIGNAL :	_____	31.7 mV
OutREC ,mV		10max

Using the **ArrowUp** button (▲) you can switch to higher maximum voltages, the **ArrowDown** (▼) button leads to smaller maximum voltages for the recorder output.

10max ▲ 100max ▲ 1000max

To accept the new settings press the **ENTER** button.

14.8 MARKERS

To modify settings concerning marker signals press the **MENU** button until the following status line is displayed:

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
UseMarker          NONE
```

The function UseMarker can be used to send a marker signal to the recorder output. Is the parameter UseMarker set to NONE no marker signal can be sent to the recorder. Pressing one of the buttons **ArrowUp** (▲) or **ArrowDown** (▼) changes the setting to **ARROW↑**.

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
UseMarker          ARROW↑
```

To accept this setting press the **ENTER** button. Is this function activated like this a marker signal can be sent to the recorder output each time the button **ArrowUp** (▲) is pressed.

The height of the marker signal is set to $\pm 1/8$ of the selected recorder signal range. When the marker signal is sent to the detector you will hear an acoustic signal.

14.9 SETTINGS FOR SERIAL PORT (RS232)

To modify settings concerning the digital communication using the serial RS232 port press the **MENU** button until the status line in the display changes to the following.

```

Temp: _____ 28.5 °C
OptBal: _____ 13.4 %%
SIGNAL: _____ 31.7 mV
ComRS232 _____ LOCK
  
```

Using this function ComRS232 you may adjust the parameters for the serial data interface. If the mode is set to "LOCK" no detector is sent to the RS232 port. Use this setting only in that case you want to record your data using an analog data system. Use the buttons **ArrowUp** (▲) and **ArrowDown** (▼) to select the frequency for data output. You can choose between two modes:

- 1 Hz One data point per second is sent to the RS232 port
- 10 Hz Ten data points per second are sent to the RS232 port

```

Temp: _____ 28.5 °C
OptBal: _____ 13.4 %%
SIGNAL: _____ 31.7 mV
ComRS232 _____ 1HZ
  
```

```

Temp: _____ 28.5 °C
OptBal: _____ 13.4 %%
SIGNAL: _____ 31.7 mV
ComRS232 _____ 10HZ
  
```

To accept your settings press the **ENTER** button. Pressing the **MENU** button will discard any new input.

14.10 MENU STRUCTURE IN NORMAL MEASURING MODE (NORMMODE)


NormMode

In this mode the user can adjust several parameter for the operation of the RI2000 refractive index detector. The NormMode is used to:

- display temperature [°C], optical balance [%] and the current detector signal [mV]
- modification and control of instrument's settings

In this mode the button **Menu** is used to switch between different sub-menus the button **Enter** is used to accept changes to parameters. The arrow buttons **◀ ▶ ▲ ▼** are used to change between several possible options.


6 seconds after switching on the RI2000 refractive index detector the logo will disappear from the display and the instrument switches to normal measuring mode.

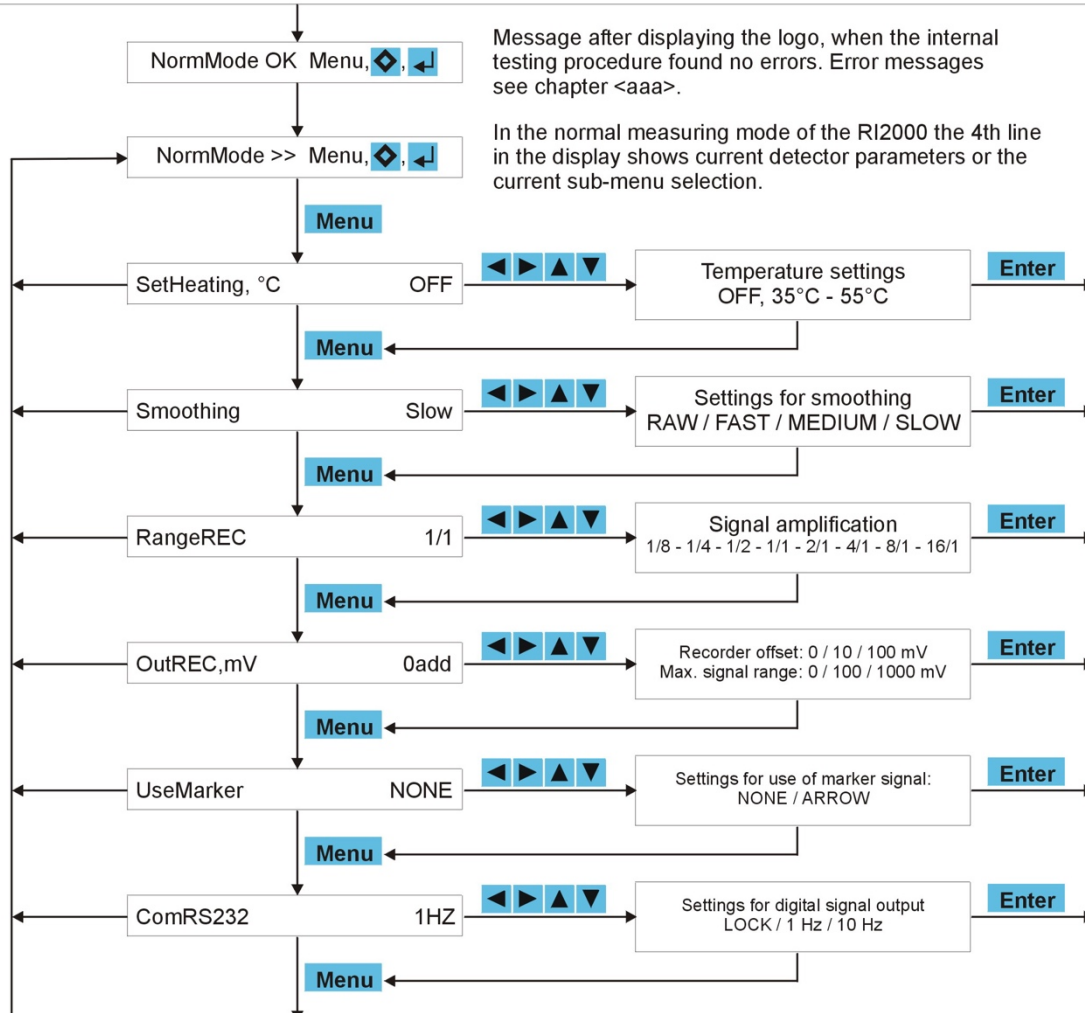

15 THE SERVICE MODE

Since version 3.x the internal firmware of the RI2000 refractive index detector contains the service mode. This mode allows the user to check more detailed system parameters for example in case of a malfunction of the instrument to get more detailed information to call the technical support.

To access service mode press the keys **ArrowLeft** (◀) and **ArrowRight** (▶) at the same time for about three seconds. The display will change from normal measuring mode (NormMode) to service mode (ServMode).

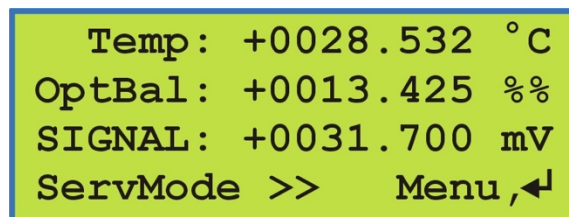

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> Menu,⏪
```

The numerical values for temperature, optical balance and the current detector signal are now displayed with three decimals. In this mode the status line is used to navigate between several sub menus.

You might choose between:

- NormMode switch back to normal measuring mode
- SwHwInfo switch to Software/Hardware information
- ViewFine switch to view fine system analysis mode
- AdjuDACs switch to adjust mode for digital analog converter
- Calibr'n switch to calibration mode

To switch back to the normal measuring mode press the **MENU** button until the NormMode statement is shown in the status line.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> NormMode
```

Now press the **ENTER** button to leave the service mode.

15.1 FIRMWARE INFORMATION

If you need detailed information about your instrument such as the firmware version number or the instrument's serial number press the **MENU** button in service mode until the menu SwHwInfo is shown in the status line.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> SwHwInfo
```

Then press the **ENTER** button. The display will show the currently running firmware version and the serial number of the instrument. The full information can be found in line two of the display.

```
S(c)hambeck SFD GmbH
RI2000 V3.2C 0401207
SIGNAL: -0031.700 mV
ServMode SwHwInfo
```

To return to the former display press the **MENU** button.

15.2 VIEWFINE MODE

To check internal settings of the refractive index detector the view fine mode can be activated. To access this mode press the **MENU** button until the sub menu ViewFine is shown in the status line.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> ViewFine
```

Now press the **ENTER** button and the display will change to the following information:

```
Check Source Unit ←
-0010 5030 2510 2520
SIGNAL: +0031.700 mV
ServMode      ViewFine
```

The values in line two represent four different voltages measured by the instrument. To get an information which voltage is shown in which column press the **ENTER** button. The display will change to the following:

```
Check Source Unit ←
Diff Summ Smpl Rfrn
SIGNAL: +0031.700 mV
ServMode      ViewFine
```

To return to the former view press the **ENTER** button again.

The displayed shortcuts represent the following voltages:

- **Diff** Difference voltage. This is the difference between the voltages resulting from measuring light intensity on the sample and the reference side of the light sensor.
- **Summ** Sum-voltage. This is the sum of the two voltages measured on the sample and the reference part of the light sensor.
- **Smpl** Sample voltage. Corresponding to the light intensity on the sample side of the light sensor.
- **Rfrn** Reference voltage. Corresponding to the light intensity on the reference side of the light sensor.

All voltages are displayed in [mV]. The value of the difference voltage can be positive or negative while the values of the other three voltages have only positive values. The sample voltage and the reference voltage can also be measured at two test points on the circuit board using a multimeter.

The sum and difference voltages are calculated according to the following formulas.

$$\text{Summ} = \text{Smpl} + \text{Rfrn}$$

$$\text{Diff} = \text{Smpl} - \text{Rfrn}$$

15.3 ADJUSTING DIGITAL – ANALOG – CONVERTERS

Press the **MENU** button until the sub menu “AdjuDACs” is shown in the display.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> AdjuDACs
```

Press the **ENTER** button to access the mode to adjust a digital to analog converter. The display will change to the following view:

```
Set both INT&REC on
fixValue: _____ mV ←
SIGNAL: +0031.700 mV
ServMode      AdjuDACs
```

This function allows the user to apply a defined signal to the recorder output and the integrator output.

■ NOTE

Make sure that the signal amplification for the recorder output is set to 1/1 before using this function. Otherwise you have to remember the amplification factor.

This voltage can be checked using a multimeter or may be used to test the used data system. By pressing the **ENTER** button you can go stepwise through a list of different voltages.

0 mV ↵ 1000 mV ↵ 1200 mV ↵ 0 mV ↵ -1000 mV ↵ -1200 mV

Pressing the **ENTER** button once will show the following information in the display:

```
Set both INT&REC on  
fixValue: 0000 mV ←  
SIGNAL: +0031.700 mV  
ServMode AdjuDACs
```

Now the signal at the integrator output and the recorder output is exactly 0 mV. Press the **ENTER** button again to switch to the next voltage. The current voltage which is applied to the signal output is shown in the display.

15.4 CALIBRATING THE RI2000

The calibration mode of the RI2000 detector is used to calibrate the signal response to a certain sample concentration. This step is needed if the instrument is used to perform quantitative analysis. The detector is calibrated when shipped to the customer. Usually a re-calibration has to be performed after modifications of the optical system only. Under normal operation conditions a re-calibration of the RI2000 is not needed.

To calibrate the signal response of the instrument you need a solution of exactly **343 mg Sucrose in 100 mL distilled water** which is used as standard. The refractive index of this solution is known and saved as an internal constant in the detectors firmware.

The first step of the calibration is purging the reference chamber and the sample chamber of the flow cell with distilled water. In the second step the sample chamber is filled with the standard solution. At this point the theoretical detector signal is 500 mV. The real detector signal is recognized by the firmware and the calibration factor to reach the theoretical 500 mV is calculated and saved as internal constant.

The correlation between the refractive index and the detector signal is:

$$1 \text{ mV} = 1 \text{ } \mu\text{RIU}$$

15.4.1 HOW TO CALIBRATE...

To activate the calibration mode press the **MENU** button until the statement “Calibr'n” is shown in the status line.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> Calibr'n
```

By pressing the **ENTER** button the calibration mode is opened. The following information will be shown in the display.

```
3-steps-Calibration:
yourTest toGet 500mV
SIGNAL: +0031.700 mV
ServMode      AdjuDACs
```

To start the calibration press the **ENTER** button. For the first step of calibration the valve will be switched into purge mode (the red purge LED will light up). The display will show the line “nowWATERwash”.

```
Step 1/3 PurgeWater
nowWATERwash, then ↵
SIGNAL: +0031.700 mV
ServMode      Calibr'n
```

At this point flush the reference chamber and the sample chamber with at least 5 mL distilled water. It is recommended to use a disposable syringe during the calibration to press the solutions into the INLET port. When the first step is finished, press the **ENTER** button. The display will change to the following view:

```
Step 1/3 PurgeWater
-0010 5030 ...WAIT !
SIGNAL: +0033.952 mV
ServMode      Calibr'n
```

The signal is recorded and checked. In the display line two the difference voltage (Diff, -0100 in this example) and the sum voltage (Summ, 5030 in this example) is shown.

If the signal is stable the display view will change to:

```
Step 2/3 InjectTest
nowTESTinto, then ↵
SIGNAL: +0471.935 mV
ServMode      Calibr'n
```

The valve is switched automatically into the measuring position. Now flush the sample chamber of the flow cell at least three times with the standard solution before pressing the **ENTER** button to proceed with step two of the calibration process. After pressing the **ENTER** button the detector registers the current signal while the display will show the different voltages as described before.

```
Step 2/3 InjectTest
0890 5060 ...WAIT !
SIGNAL: +0471.935 mV
ServMode      Calibr'n
```

If the signal is stable the calibration factor is calculated and saved. The success of the calibration process is confirmed by a displayed message.

```
Step 3/3 SaveSession  
Calibration was OK ↵  
SIGNAL: +0492.980 mV  
ServMode      Calibr'n
```

Press the **ENTER** button to finish the calibration process. The new calibration factor is shown in the upper right corner of the display. The four measured and calculated voltages are shown in the display according to the ViewFine mode (see [chapter 15.2](#)). The calibration is now finished.

```
set new SPAN  02803  
-0010 5030 0890 5060  
SIGNAL: +0492.980 mV  
ServMode      Calibr'n
```

In case of a problem during the calibration, an error message will be displayed.

```
Step 3/3 SaveSession  
Calibration failed  
SIGNAL: +0492.980 mV  
ServMode      Calibr'n
```

In this case the RI2000 uses the old calibration factor for further measurements. Repeat the process of calibration as describe before.

15.5 SERVMode OVERVIEW

ServMode

This mode allows the user to get detailed information about the current status of the RI2000. If the calibration mode is not entered no changes can be applied to the instrument.


Features of the service mode:

- more detailed display of temperature (°C), optical balance (%%) and detector signal (mV)
- system information (such as firmware version and serial number of the instrument)
- display of optical bench parameters (reference / sample voltages, sum- and difference voltage)
- calibration of the detector response (3-step-calibration)

In the service mode only the buttons **Menu** and **Enter** are active to select the different sub-menus and to accept new settings.

To discard changes and to leave a sub-menu press the **Menu** button.

By pressing the buttons ArrowLeft ◀ and ArrowRight ▶ for three seconds at the same time the instrument will switch into **ServMode**

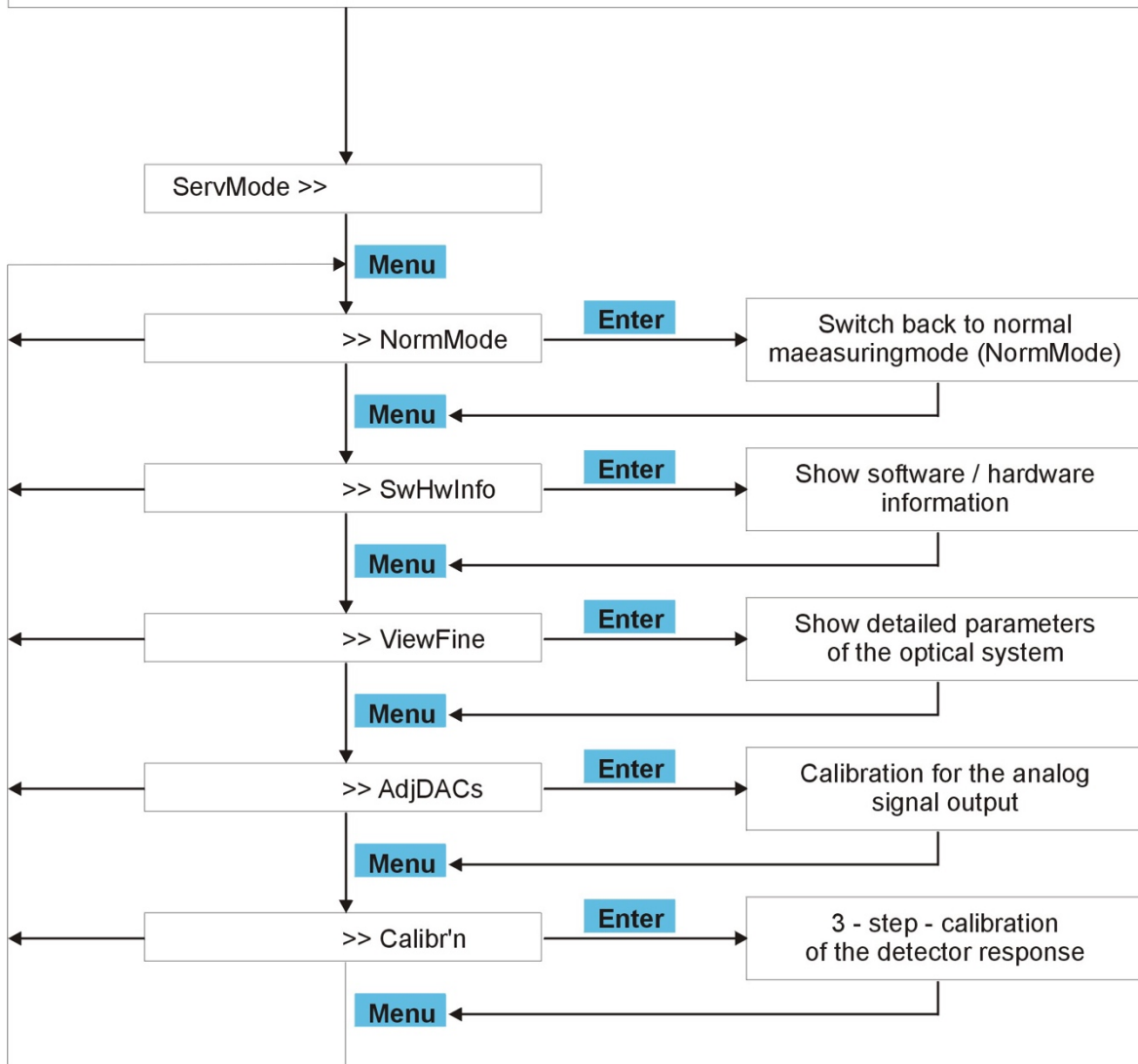


15.6 SERIAL COMMUNICATION

REFRACTIVE INDEX DETECTOR RI2000 - SERIAL RS232 INTERFACE - FIRMWARE V. 3.2

The RS232 port is used for the communication of detector and computer.

The following options can be used to acquire data and control the RI2000 refractive index detector:

- AQUIRE2000 - software solution for data acquisition and instrument control
- special designed software solutions depending on the customer's needs
- simple software solutions like Windows HyperTerminal or special macros for Microsoft Excel

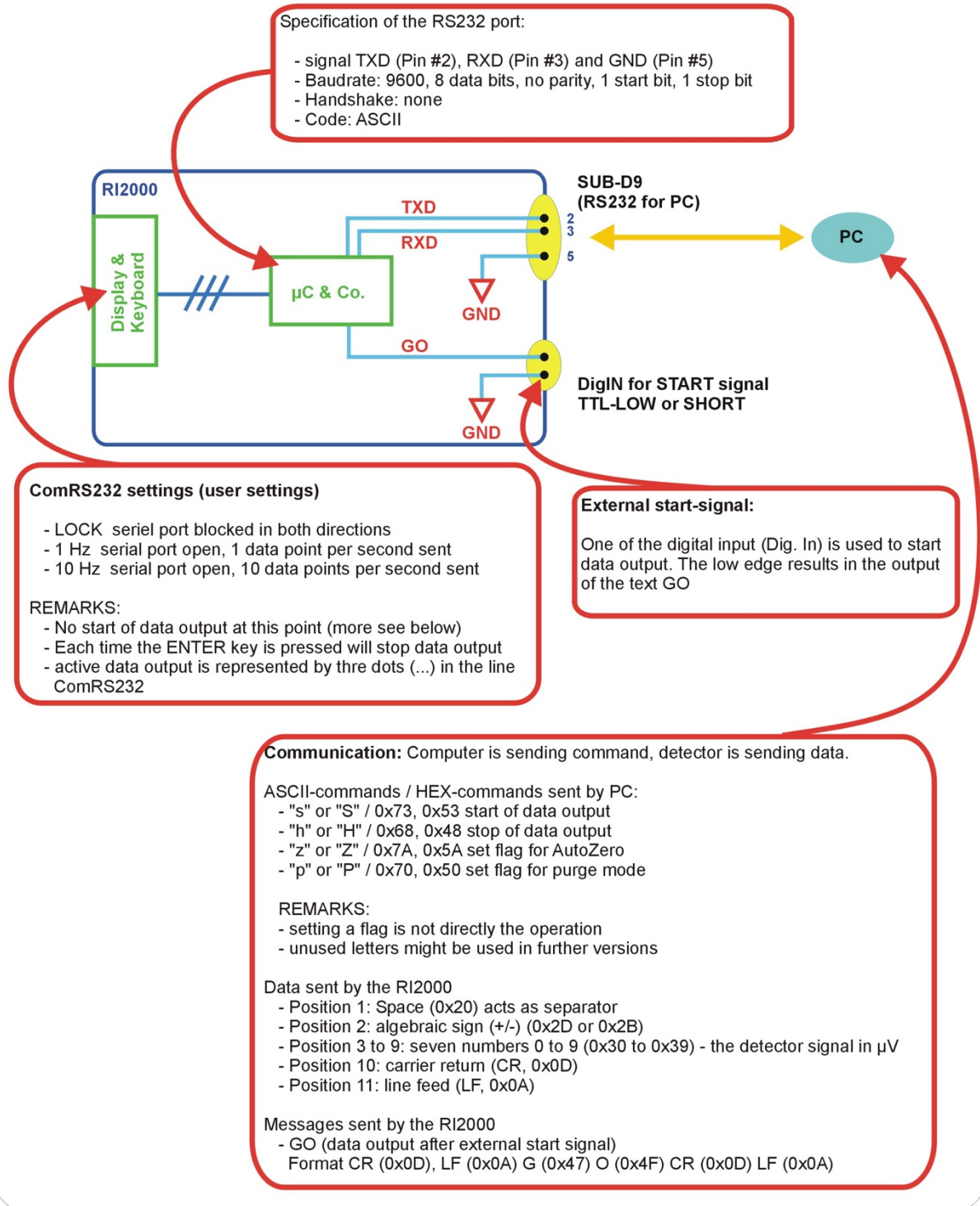

16 MAINTENANCE OF THE RI2000

16.1 IDENTIFYING HARD- AND SOFTWARE

If you need a detailed information about your instrument such as the firmware version number or the instrument's serial number you can get this information in the service mode of the detector. Press the buttons **ArrowLeft** (◀) and **ArrowRight** (▶) at the same time for approx. three seconds. The display will switch into service mode.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> Menu,⏪
```

Press the **MENU** button until the menu SwHwInfo is shown in the status line.

```
Temp: +0028.532 °C
OptBal: +0013.425 %%
SIGNAL: +0031.700 mV
ServMode >> SwHwInfo
```

Then press the **ENTER** button. The display will show the currently running firmware version and the serial number of the instrument. The full information can be found in line two of the display.

```
S(c)hambeck SFD GmbH
RI2000 V3.2C 0401207
SIGNAL: -0031.700 mV
ServMode SwHwInfo
```

To return to the former display press the **MENU** button. To exit the service mode select the menu NormMode by pressing the **MENU** button and press **ENTER**.

16.2 ERROR MESSAGES

In case of some malfunction the RI2000 refractive index detector may display an error message. The most error will appear directly after switching on the instrument.

16.2.1 MALFUNCTION OF THE TEMPERATURE SENSOR

There are two possible error messages indicating problems with the temperature sensor. A more detailed description of the temperature control of the RI2000 can be found in [chapter 16.8](#).

The message “noTS” (no temperature sensor) will appear for that reason that the temperature sensor is not connected to the circuit board or not working.

```
Temp: _____ noTS °C
OptBal: _____ 13.4 %%
SIGNAL: _____ 31.7 mV
NormMode >> Menu, ⬠, ⬅
```

The message “!” is shown if the temperature of the optical bench is outside the range between 9 °C and 64 °C. Beside the displayed message an acoustic signal will occur.

```
! Temp: _____ 8.2 °C
OptBal: _____ 13.4 %%
SIGNAL: _____ 31.7 mV
NormMode >> Menu, ⬠, ⬅
```

16.2.2 MALFUNCTION OF THE ANALOG DIGITAL CONVERTER

This error messages only occur directly after switching on the RI2000 refractive index detector.

```
Temp: _____ 28.5 °C
OptBal: _____ noAD %%
SIGNAL: _____ 31.7 mV
NormMode >> Menu, ⬠, ⬅
```

If this message “noAD” shows up the analog digital converter did not connect to the micro controller. Switch of the instrument, wait for ten seconds and try to start again. If the problem remains contact the technical support.

```
Temp: _____ noTS °C
OptBal: _____ noAD %%
SIGNAL: _____ mV
ServMode >> NormMode
```

After switching to service mode you will also notice the error message “noTS”. As the micro controller also controls the temperature of the optical bench this error message shows up but does not mean that there is a problem with the temperature sensor.

16.2.3 INTENSITY ALARM

The intensity alarm message can show up directly after initialization of the instrument or during operation. There are two different ways of displaying this message.

```
Temp: _____ 28.5 °C
OptBal: _____ 13.4 %%
SIGNAL: intensity mV
NormMode >> Menu, ⬠, ⬅
```

```
Temp: _____ 28.5 °C
OptBal: _____ 13.4 %%
SIGNAL: INTENSITY mV
NormMode >> Menu, ⬠, ⬅
```

It is normal that you get no beep from the RI2000 refractive index detector when you press the **AutoZero** or **PURGE** button while the 'INTENSITY' or 'intensity' error message is shown. In this case it makes no sense to press the **AutoZero** or **Polarity** buttons when there is no signal value shown in the display. The functions **AutoZero** and **Polarity** are only active if there is a detector signal.

Is the message "intensity" displayed the voltage U_{Summ} is below 1.7 V. The reason is low light intensity reaching the light sensor. To resolve this problem, try the following:

Purge the reference chamber and the sample chamber of the flow cell and make sure that both chambers are filled with mobile phase

- Check the light source (see chapter 16.4)
- Make sure the flow cell is clean
- Adjust the alignment of the optical system (see chapter 16.3)
- Check the lamp current, it might be too low.

In case of the message "INTENSITY" (all capital letters) the problem is that the voltage U_{Summ} is beyond 7.5 V. This might be caused by several facts:

- The lamp current is too high, check lamp current
- The cover of the optical system is open and light from outside reaches the light sensor

16.2.4 SIGNAL OUT OF RANGE

The error message of the type "OVER_ADxx" shows up if the analog digital converter gets a signal which is too high or too low.

```

Temp:      _____ 28.5 °C
OptBal:    _____ 13.4 %%
SIGNAL:    OVER_AD1H mV
NormMode >> Menu, ⬠, ⬅

```

```

Temp:      _____ 28.5 °C
OptBal:    _____ 13.4 %%
SIGNAL:    OVER_AD1L mV
NormMode >> Menu, ⬠, ⬅

```

```

Temp:      _____ 28.5 °C
OptBal:    _____ 13.4 %%
SIGNAL:    OVER_AD2H mV
NormMode >> Menu, ⬠, ⬅

```

```

Temp:      _____ 28.5 °C
OptBal:    _____ 13.4 %%
SIGNAL:    OVER_AD2L mV
NormMode >> Menu, ⬠, ⬅

```

To localize the problem the “xx” in the error message is filled with a two letter code according to the following table

xx=1H	sample voltage beyond 5 V
xx=2H	reference voltage beyond 5 V
xx=1L	sample voltage below 0 V
xx=2L	reference voltage below 0 V


The sample voltage corresponds to the light intensity reaching the light sensor on the sample side, as the reference voltage corresponds to the light intensity reaching the light sensor on the reference side.

When this error message shows up, try the following to solve the problem:

- Purge the flow cell to ensure that both chambers are filled with mobile phase (without bubbles)
- Check the optical balance ([see chapter 16.3](#))

16.2.5 RESET CONFIGURATION

This type of error might only occur when the instrument is switched on. If the error message “RC” is shown in the display the controller was not able to load the internal configuration during the initialization.

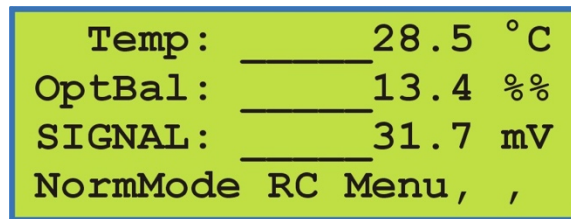

```
Temp : _____ 28.5 °C
OptBal : _____ 13.4 %%
SIGNAL : _____ 31.7 mV
NormMode RC Menu , ,
```

In this case for example the calibration factor is not loaded. The instrument's setting will be set automatically to factory default settings. If you are performing quantitative analysis you have to re-calibrate your detector before you can proceed.

16.3 ADJUSTING THE OPTICAL SYSTEM

By adjusting the optical system of the RI2000 refractive index detector it is possible to influence the light intensity reaching the two photodiodes of the light sensor.

To perform the adjustment the sample chamber and the reference chamber of the flow cell have to be purged with distilled water. When both chambers contain the same liquid theoretically the same light intensity should reach the sample side and reference side of the light sensor. In this case the optical balance should be zero.

The adjustment of the optical system is performed manually. During this procedure the position of the light beam on the photo sensor is changed by adjusting the mirror.

During the procedure you can follow the changes by observing the change of the detector signal or more detailed by observing the difference and sum voltage.

For the calculation of the different voltages, the detector signal and the optical balance the following formulas are used:

$$Diff = Smpl - Rfrn$$

$$Summ = Smpl + Rfrn$$

$$OB = \frac{Diff}{Summ} \cdot 1000 (\text{‰})$$

$$Signal = POL \cdot (OB \cdot C_0 - AZ)$$

Where Signal is the calculated detector signal, POL the selected polarity, C_0 is the calibration factor resulting from the 3-step-calibration and AZ the 'real' signal before pressing the AutoZero button. The pre-set calibration factor is 2724. This factor is used if no 3-step-calibration of the instrument was performed.

To adjust the optical system follow the steps listed below.

- Switch to purge mode, make sure the red purge LED is on.
- Flush both chambers of the flow cell with clean mobile phase for several minutes.

- Switch back to measuring mode, make sure the red purge LED is off.
- Remove the black caps at the back side of the instrument
- Put the hex-wrench into one of the opened hole

■ NOTE

There are two adjustment screws for the mirror. Keep in mind that you might have to turn the other screw. Anyway you should turn the screws not more than $\frac{1}{2}$ or $\frac{3}{4}$ turns.

- Adjust the screws until the optical balance reaches a value of $0 \pm 10 \%$.

Now the adjustment for the used solvent is finished.

16.4 LAMP EXCHANGE AND ADJUSTMENT OF THE LAMP

It might be possible that you need to re-adjust the light source due to some changes in the used light bulb due to the transport or you need to replace a burned out lamp. The adjustment of the light source is only possible when the housing is opened. To prepare for adjustment do the following steps:


- Switch off the instrument
- Unplug the mains cable to prevent electrical shock when opening the housing

■ NOTE


Some operations during the adjustment have to be done when the instrument is in operation. So you need to re-connect the mains cable and to switch on the RI2000 refractive index detector. When the instrument is operated with opened housing make sure that you do NOT touch any electrical component!




- Open the outer housing of the optical bench (chrome coloured)
- Open the inner black housing of the optical bench. You will see the optical system of the RI2000 refractive index detector. The components are shown below:

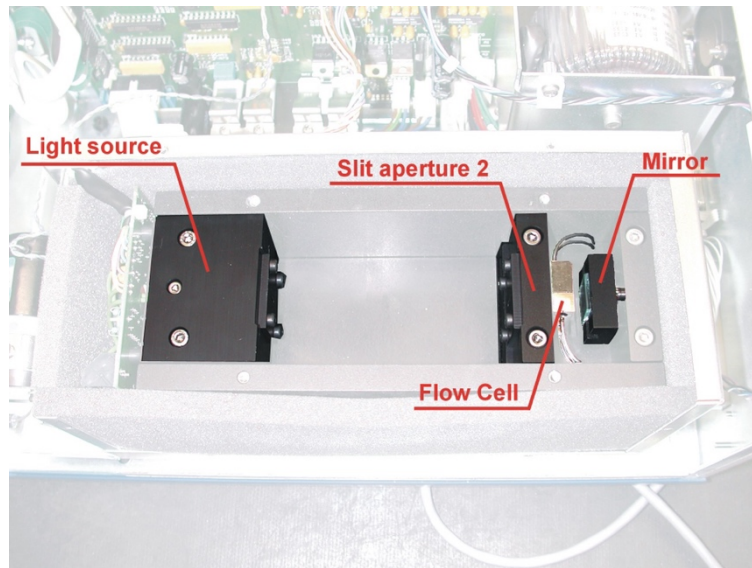

- Loosen the holding screw of the light source and remove the lamp.

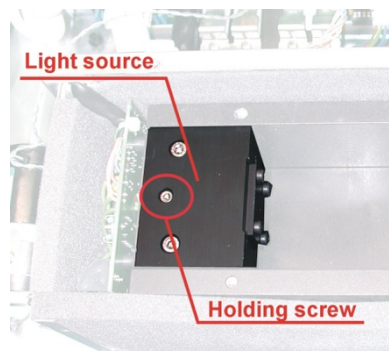

- Loosen the screw connector of the lamp's power supply at the pre-amplifier board (labelled: LAMPE)

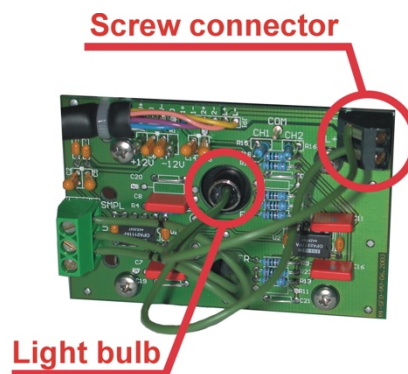

- Place a new lamp in the holder and tighten the holding screw

- Connect the power wires of the lamp with the screw connector on the pre-amplifier board.
- Connect the mains cable to the detector and switch on the instrument

■ **NOTE**

The following steps need to be performed on the running instrument. There is the danger of electrical shock when touching electrical components.

- Use a voltmeter to check the lamp voltage at the test points on the circuit board. The label of this test point depends on the version of this board. If the serial number of your instrument is 0312xxx the points are labelled 'Lampe' and 'GND', if the serial number of your instrument is 0401xxx the points are labelled LampeULmp und 'GND'. To find your serial number see [chapter 15.1](#). The lamp voltage should be 3.3 ± 0.3 V. If needed you can adjust this voltage by turning the potentiometer R19 at the main circuit board.


In the following two steps the light bulb is positioned for optimal operation.

Step 1:


Place the lamp in the holder resulting in a sharp picture of the illuminated area at the slit aperture 2. Usually the metal cover of the light bulb ends with the holder, sometimes the metal cover stands about 1 or 2 mm out of the metal block.

Step 2:

Turn the light bulb that the picture of the filament is parallel to the slit of the second slit aperture.

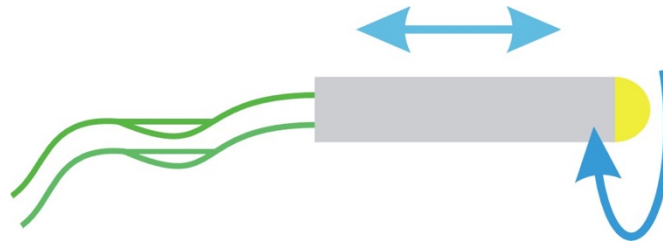


Make sure that only 2/3 of the slit are illuminated by the light source. This is illustrated in the following figure.

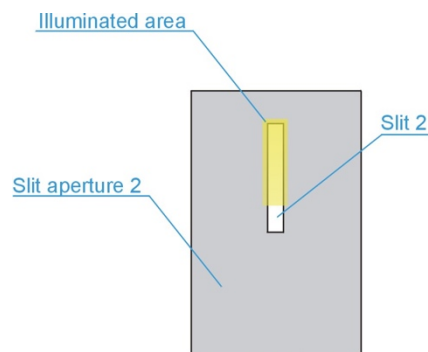

■ NOTE


The projected picture at the light sensor should be a sharp one. The full height of the photo diodes has to be illuminated by the light beam. This is shown in the following figure.

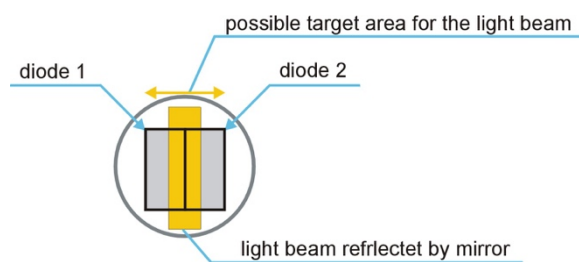

16.5 CHECKING SUM- AND DIFFERENCE VOLTAGES

This chapter describes the check of the sum and difference voltages. These voltages result from different light intensities reaching the sample side and the reference side of the light sensor. To perform the check follow the steps listed below:

- Make sure the instrument is switched on and both chambers of the flow cell are filled with distilled water (without bubbles). In this case the optical balance should be $0 \pm 20,0 \%$. To ensure a stable temperature the instrument should be switched on for several hours before performing this test.
- Press the **PURGE** button. Make sure the purge mode is activated, the red purge LED should light up.
- Use a disposable syringe to press approx. 5 mL distilled water through the flow cell
- The display will show the current detector signal. If no signal is displayed please check chapter **XXX**.
- Activate the service mode of the RI2000 by pressing the buttons **ArrowLeft** (◀) and **ArrowRight** (▶) at the same time for approx. three seconds. The display will change to the following

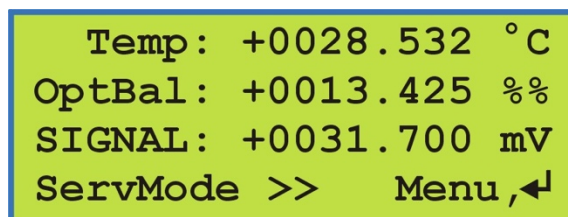

The image shows a screenshot of the instrument's LCD display. The display is green and shows the following text: Temp: +0028.532 °C, OptBal: +0013.425 %, SIGNAL: +0031.700 mV, and ServMode >> Menu, ←.

- Press the **MENU** button until the status line in the display changes to ViewFine. Now press the **ENTER** button. The display will show four values.

```
Check Source Unit ←  
-0010 5030 2510 2520  
SIGNAL: +0031.700 mV  
ServMode      ViewFine
```

To find out which value is displayed in which column press the **ENTER** button.
The display will change to the following view:

```
Check Source Unit ←  
Diff Summ Smpl Rfrn  
SIGNAL: +0031.700 mV  
ServMode      ViewFine
```

To go back to the former view press the **ENTER** button again.

Diff represents the difference voltage (in mV) and corresponds to the difference of intensity reaching the sample side and the reference side of the light sensor.

Summ represents the sum voltage (in mV) and corresponds to the light intensity reaching both sides of the light sensor.

Smpl represents the sample voltage which corresponds to the light intensity reaching the sample side of the light sensor.

Rfrn represents the reference voltage which corresponds to the light intensity reaching the reference side of the light sensor.

The voltages are calculated according to the following formulas

$$\text{Diff} = \text{Smpl} - \text{Rfrn}$$

$$\text{Summ} = \text{Smpl} + \text{Rfrn}$$

- Use a multimeter to check the lamp voltage at the test point on the circuit board. For this you need the serial number of your instrument ([see chapter 15.1](#)) to locate the test points. If your serial number is 0312xxx the points are labelled “Lampe” and “GND”. In case of a serial number of 0401xxx the points are labelled “LampeULmp” and “GND”. In both cases the lamp voltage should be set 3.3 ± 0.3 V. On pages 72 and 73 you can find a schematic view of the main circuit board with the marked positions of the test points.
- If the sum voltage is only a few millivolt the light bulb does not work. Check the power cables of the light source to make sure that they are connected properly. If needed check the light source ([see chapter 16.4](#)).
- The sum voltage should be in the region of 5000 ± 500 mV. If this is not the case adjust the lamp voltage by turning the potentiometer R19 until the sum voltage is in the region of 5000 ± 500 mV.
- Due to the adjustment of the light intensity it is necessary to check the lamp voltage again. If this voltage is not in the region von 3.3 ± 0.3 V the light source has to be re-adjusted.
- The optical system is adjusted properly if the sum voltage is 5000 ± 500 mV and the lamp voltage 3.3 ± 0.3 V.

16.6 CHECK AND REPLACEMENT OF THE VALVE

If the valve which is used to switch the path of the mobile phase in the normal measuring mode and purge mode needs to be replaced follow the following steps:

- Switch of the RI2000 refractive index detector and unplug the mains cable

■ **NOTE**

To prevent electrical shock make sure that the power cable is disconnected before opening the housing of the instrument.

- Remove the stainless steel capillaries from the valve body.

■ **NOTE**

Remind the position of the capillaries for proper re-connection.

- Loosen the control cable of the valve at the main circuit board.
- Loosen the two screws holding the valve body.
- Remove the valve body and replace with the new one.
- Tighten the holding screws and reconnect the capillaries in the correct positions.
- Connect the signal cable to the contact labelled "Valve" on the main circuit board.

- Start you chromatography pump to flush mobile phase through your RI2000 detector.
- Check the capillary connections for any leakage.
- Close the instrument's housing.

16.7 CHECKING AND CLEANING THE FLOW CELL

In some cases it might be necessary to clean the flow cell inside the RI2000 refractive index detector. You should try to wash away possible contaminations by purging both chambers of the flow cell with fresh mobile phase for a longer time.

If you are not sure if the problems are caused by a contaminated flow cell or caused by other problem contact the technical service of your local distributor or SRI Instruments Europe GmbH directly.

Possible reasons to open the optical bench for checking the flow cell:

- Drop of the sum voltage below 4500 mV. Make sure that your sample chambers are washed and filled with distilled water (without bubbles) before checking the sum voltage.
- Noisy baseline (may be caused by small particles in the flow cell)
- Constant drift of the baseline

To check the flow cell inside your detector follow the steps listed below:

- Switch off the instrument and unplug the power cable. Open the detector's housing after loosening the screws.

NOTE

When performing service work inside the instrument's housing make sure that the power cable is disconnected.

- Open the outer housing of the optical bench after loosening the screws
- Remove the heat insulation on top of the optical bench

- Open the inner (black) housing after loosening the screws
- Loosen the two hex-screws holding the slit aperture 2 and remove the aperture carefully
- Use a flash light to check the flow cell for contaminants or damages. Make sure that no bubbles are inside the flow cell.
- If there are contaminants in the flow cell purge it with a suitable solvent. If you are working with aqueous systems distilled water is convenient. If you use organic substances try solvents like acetone, tetrahydrofuran or chloroform. After purging the cell with solvent to clean for a longer time you should purge it with your mobile phase and check again if the contaminants is removed.

■ NOTE

In the case of dangerous solvents make sure that they are disposed correctly.

■ NOTE

In case of aqueous solvents it is possible that algae grow inside your chromatographic system. For that reason it is not recommended to store the detector for a longer time with aqueous solvents inside. For long time storage it is recommended to purge the flow cell with ethanol followed by air.

If possible we recommend the addition of a small amount of organic solvent (such as isopropanol or methanol) to your mobile phase to prevent algae growth.

If it is not possible to remove the contaminant from your flow cell it might be necessary to replace it with a new flow cell. Please contact your local distributor or SRI Instruments Europe GmbH directly for assistance.

If the flow cell is cleaned set the slit aperture back in it's position and fix it
Put the lid on the inner housing and tighten the screws

Place the heat insulation back on top of the inner housing

Close the outer housing and tighten the screws

Perform a measurement to check if the problem is solved now. If you still have problems using your RI2000 refractive index detector please contact you local distributor or SRI Instruments Europe GmbH directly for assistance.

16.8 THE HEATING CIRCUIT OF THE RI2000

In the following chapter the complete specification of the RI2000 heating circuit is put together.

- In the normal measuring mode (NormMode) the temperature is displayed with one digit (e. g. 27.1 °C). If you want a more detailed information you can switch to the service mode (ServMode) where the temperature is displayed with three digits (e.g. 27.094 °C).
- The temperature sensor will be detected automatically when the instrument is switched on.
- If the heating is switched off the optical bench is operated at a temperature of about 6 °C above ambient temperature due to heat irradiated of electronic compounds.
- In the interval from 35 °C up to 55 °C the temperature can be selected in 1 °C steps. The temperature difference of 6 °C has to be considered. If the heating is activated a asterisk (*) is shown in the upper left corner of the display.
- The heat cartridge used in the RI2000 refractive index detector has an internal resistance of 750 Ω .
- The temperature sensor is controlled by the micro controller. Heating is deactivated and cannot be activated for one of the following reasons:
 - The temperature sensor does not work or is not connected to the main circuit board. In this case the error message “noTS” will be displayed.
 - The current temperature of the optical bench is below 9 °C or beyond 65 °C. In this case the error message “!” will show up and an acoustic signal occurs.
 - If the heating is not switched off properly by the firmware a thermal fuse will switch off the heating at a temperature of 72 °C.

17 VOLTAGE TEST POINTS

The following list contains different test points on the main circuit board of the RI2000 refractive index detector.

All voltages listed below are measured in relation to Ground GND or AGND. Both are connected to the instrument's housing through the holding screws on the circuit board. Check the serial number of your instrument ([see chapter 15.1](#)) to locate the test points on the board. The figure on page 74 shows the board used in instruments of the 0312xxx series. The figure on page 75 shows the board used in instruments of the 0401xxx series.

GND	digital ground
AGND	analog ground

Uncontrolled voltages

U _d *	10.8 V	for RI2000 with serial number up to SN 0212220
	11.3 V	for RI2000 with serial number from SN 0212221

Stabilized voltages

V _{cc}	controlled voltage for digital circuits
U _{Lmp}	controlled lamp voltage
U _{ref} or TP14	reference voltage for the A/D-converter

Test point	Voltage	Adjustment
V _{cc}	5.0 V	-
U _{Lmp}	2.5 to 5.1 V	Potentiometer R19
+12 V	+12 V	-
-12 V	-12 V	-
+5 V	+5 V	-
-5 V	-5 V	-
U _{ref} or TP14	2.5 V	-

Signal voltages

The voltages you find at this test point should correspond to the voltages shown in the detector's display in ViewFine mode.


Test point	Voltage
ch1	corresponding to Smpl voltage of the optical bench
ch2	corresponding to Rfrn voltage of the optical bench
INT	voltage of the analog integrator output
REC	voltage of the analog recorder output

18 SPARE PARTS


Article No.	Description
RI2000-001	Flow cell for the RI2000 refractive index detector
RI2000-002	Valve for the RI2000 refractive index detector
RI2000-012	Tungsten lamp for the RI2000 refractive index detector
RI2000-023	Fan for the RI2000 refractive index detector
RI2000-024	Set of sealings for the RI2000 refractive index detector
90-ZBU1	Outlet port screws (0.75 mm I.D.) for RI2000
90-ZBU1C	Inlet port screws (0.25 mm I.D.) for RI2000

Spare parts you can order from your local distributor or from SRI Instruments Europe GmbH directly.


19 TEST POINTS FOR LAMP VOLTAGE (FOR SERIAL NUMBER 0312XXX)

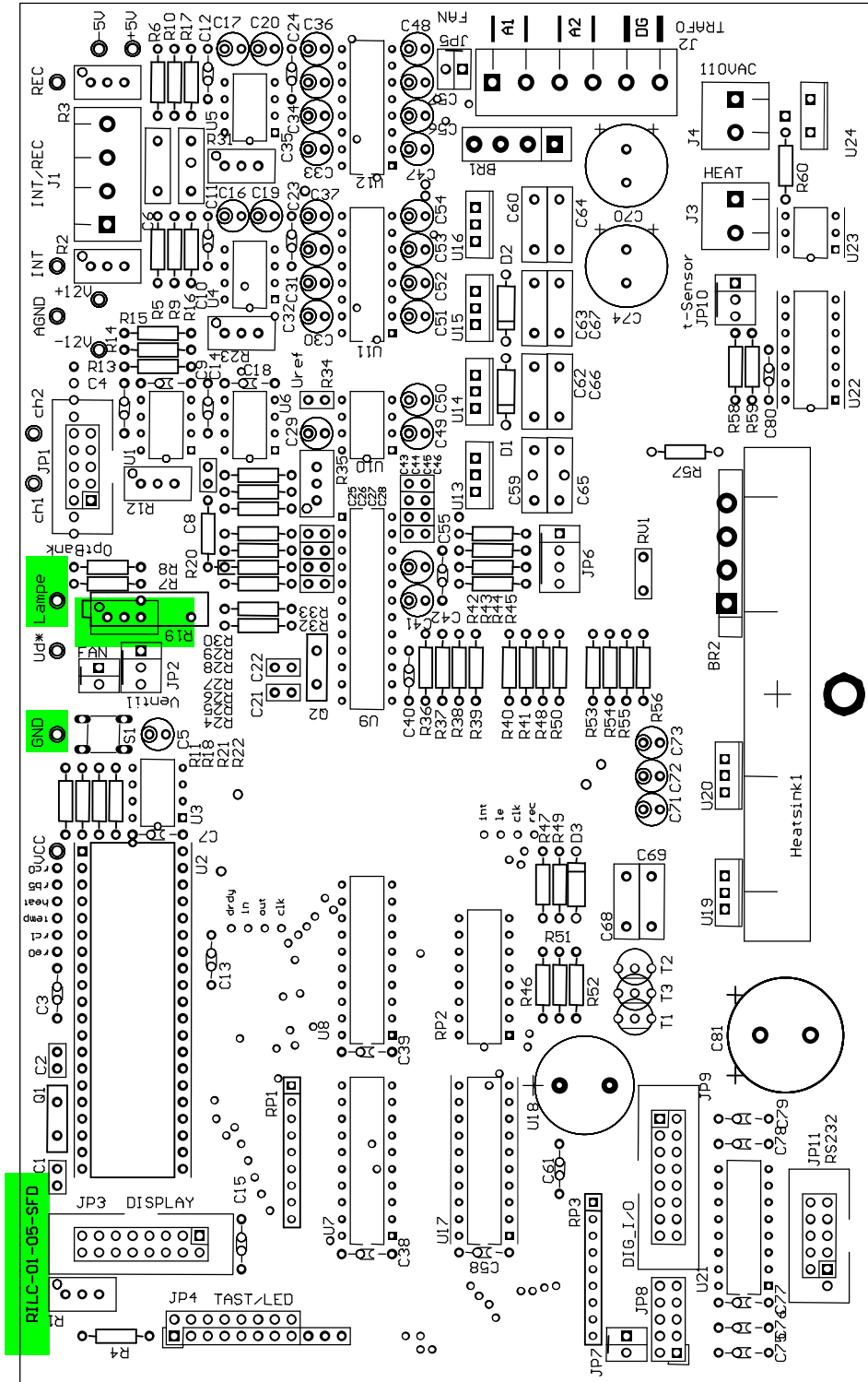

20 TEST POINTS FOR LAMP VOLTAGE (FOR SERIAL NUMBER 0401XXX)

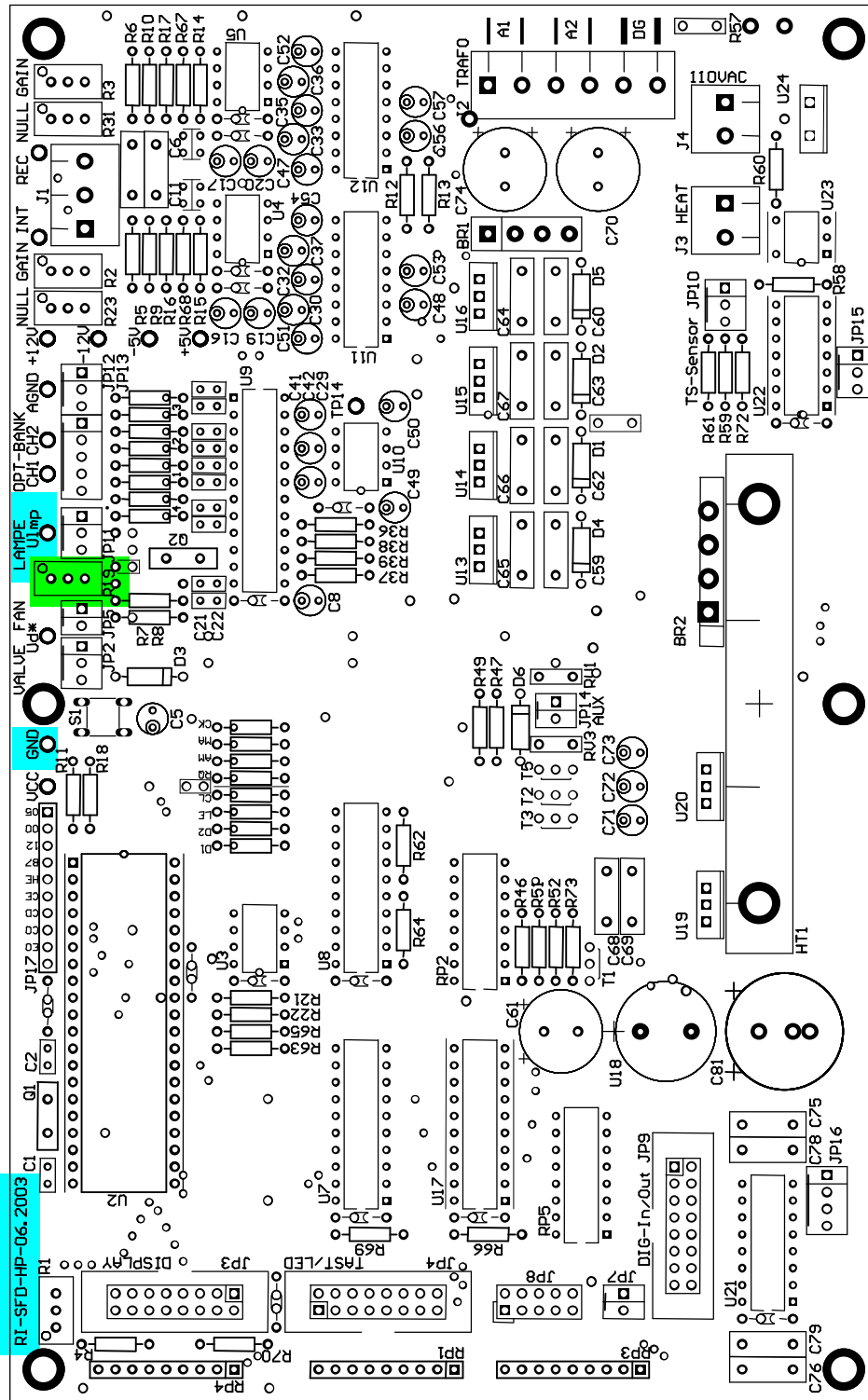

21 VOLTAGE TEST POINTS (OVERVIEW) FOR 0312XXX SERIES

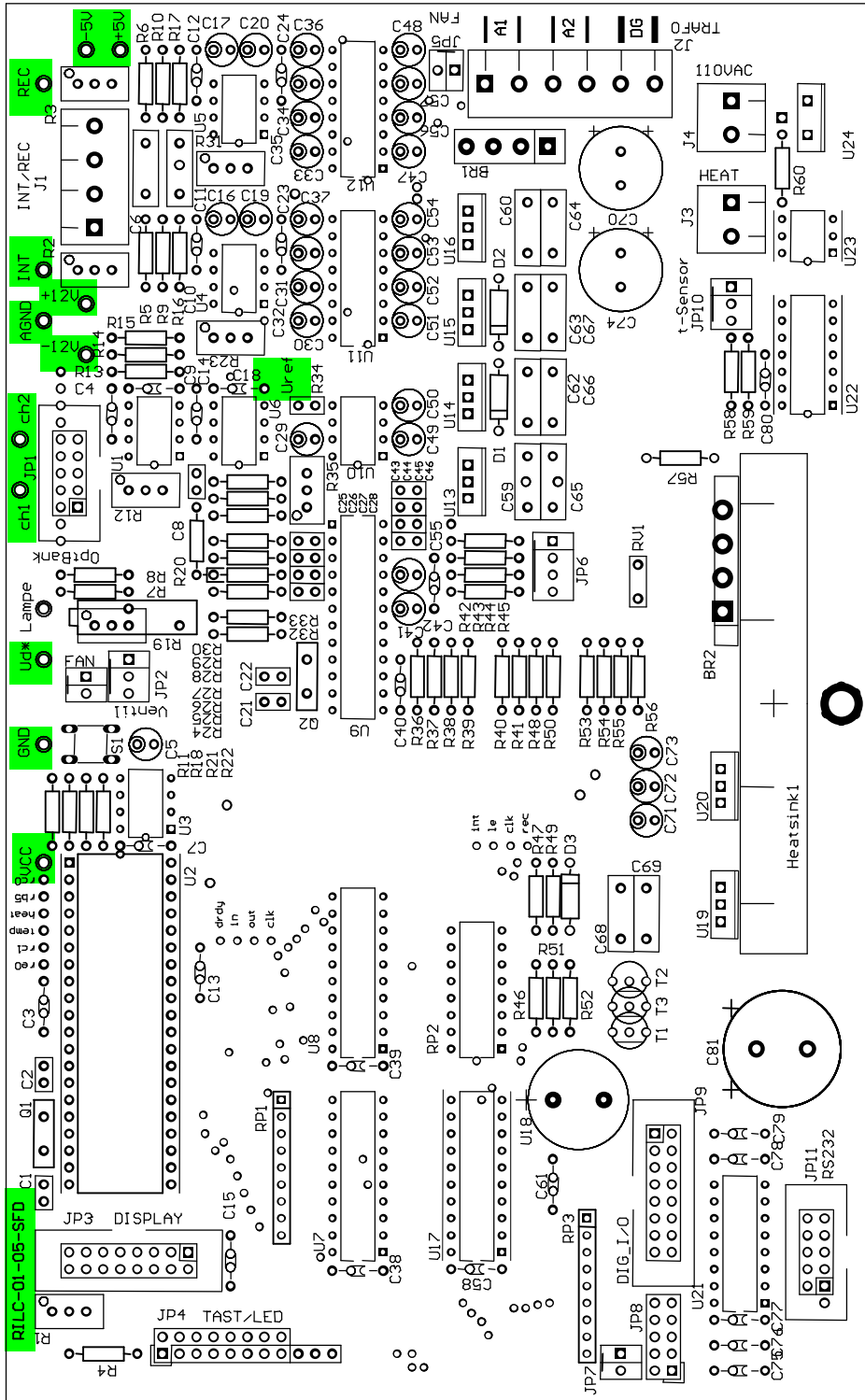

22 VOLTAGE TEST POINTS (OVERVIEW) FOR 0401XXX SERIES

